

To the users of Cios Alpha VA30 systems

Name: Andreas Fischer
Department: SHS AT SU M&S-M
E-Mail: andreas.anf.fischer@siemens-healthineers.com
Date: February 27, 2019

Important customer safety notice regarding corrective field action:

AX014/19/S

Important customer safety notice AX014/19/S regarding Cios Alpha VA30

Dear customer,

Your Cios Alpha is a high-quality mobile C-arm, which offers you excellent image quality with high generator performance and extensive protection against overheating at high power levels. Nevertheless, we would like to inform you that with the Cios Alpha, permanent operation at very high tube outputs results in increased wear of the radiation generation components.

What is the problem to be corrected and when does it occur?

The organ programs pre-installed ex works enable safe operation within the specified performance limits.

However, it is possible to overwrite the organ programs with higher performance parameters.

How does the problem affect system operation and what is the potential risk?

Permanent operation at very high tube output results in increased attrition of the radiation generation components.

This is the case, for example, when tube voltages of 125 kV are used over a longer period of time, e.g. more than 30 minutes. There is a risk of premature failure of the monoblock and thus a loss of the imaging X-rays.

What measures are being taken?

We therefore recommend avoiding overwriting the organ programs, as this can lead to more frequent operation with the performance parameters mentioned above. Please read the operating instructions in general and in this context especially the new Addendum which we have now delivered to you. Please forward this information to your employees or organizations that are affected.

How was the problem detected?

The problem was detected during our field observation.

How is the corrective measure implemented?

This letter and the addendum to the user manual will be distributed as an AX014/19/S update to Cios Alpha VA30 to affected customers.

What are the risks for patients who have been previously examined/ treated with this system?

In this case, we do not consider a follow-up examination of patients to be necessary. The problem at hand here is a possible hardware error without effect on a previous diagnosis and treatment of patients.


We would like to thank you for your cooperation in dealing with this safety notice and ask you to immediately pass this information on to all employees in your institution who need to know about this problem and to instruct them accordingly. Please also forward this safety notice to other entities that may also be affected by this measure.

If the device has been sold and is therefore no longer in your possession, we would like to ask you to forward this safety note to the current owner. If possible, please inform us of the identity of the current owner.


We thank you for your attention and cooperation and wish you continued success with your Cios Alpha.

Best regards,

SIEMENS Healthcare GmbH
Business Area Advanced Therapies

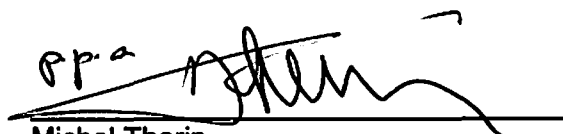

Michel Therin
Head of Business Area AT

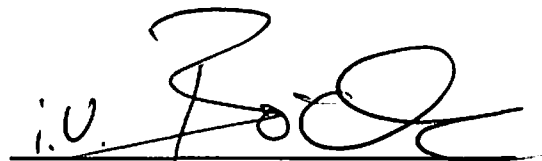

Johann Böck
Safety Officer Medical Devices AT

Important additional information

Please observe the following additional information for the Cios Alpha system.

Safety

Overload protection



Caution

Under extreme operating conditions (continuous operation close to the maximum tube voltage), the X-ray tube can become very hot. This can shorten the service life of the system components.

Risk of burns and prolongation of the operation!

- ◆ Avoid touching the X-ray tube housing.
 - ◆ Change your parameter settings, e.g. by reducing the dose level.
 - ◆ Please make sure that the system has cooled down before performing critical interventional applications.
 - ◆ We also recommend keeping a second system ready for continuous radiation lasting longer than 30 minutes.
-

Important additional information

