

Opfikon, 20. Dezember 2023

Dringender Sicherheitshinweis: FA-2023-055: Q-stress 6 und Xscribe 6

Sehr geehrte Damen und Herren,

Baxter gibt eine Korrektur für die Herzbelastungstestsysteme Q-Stress und XScribe (Q-Stress Version 6 oder höher und XScribe Version 6 oder höher) aufgrund einer möglichen Änderung der QRS-Amplitude aus, die in Elektrodiagramm-(EKG)-Messwerten festgestellt wurde. Der Quellkonsistentfilter (SCF) ist aktiviert. Der SCF ist eine konfigurierbare Funktion bei Mehrkanalaufzeichnungen, die darauf abzielt, Rauschen zu reduzieren und inkonsistente Signale zu entfernen, ohne das mit Stresstests verbundene EKG-Signal zu verzerren.

Baxter führte interne Tests anhand der EKG-Daten des Patienten mithilfe des SCF-Algorithmus durch, die eine durchschnittliche QRS-Amplitudenänderung von 5,4% Reduzierung vor dem Training und 7,1% Reduzierung während der Spitzenherzfrequenz zeigten (siehe Abbildung 1). Die beobachteten maximalen QRS-Amplitudenänderungen waren eine Reduzierung um 39,7% während der Spitzenherzfrequenz.

Wenn der SCF aktiviert ist, können QRS-Amplitudenänderungen in der angezeigten EKG-Kurve auf der Echtzeitanzeige auf dem Bildschirm, in Life-EKG-Ausdrucken und in Abschlussberichten beobachtet werden. Die Anzeige des durchschnittlichen Herzschlags und alle Berechnungen (z.B. Herzfrequenz, ST-Level, ST-Steigung) basieren nicht auf SCF-gefilterten Daten; Daher sind sie nicht betroffen. Durch Algorithmen erkannte Ereignisse (z.B. PVC, VRUN) basieren ebenfalls nicht auf SCF-gefilterten Daten und sind daher nicht betroffen.

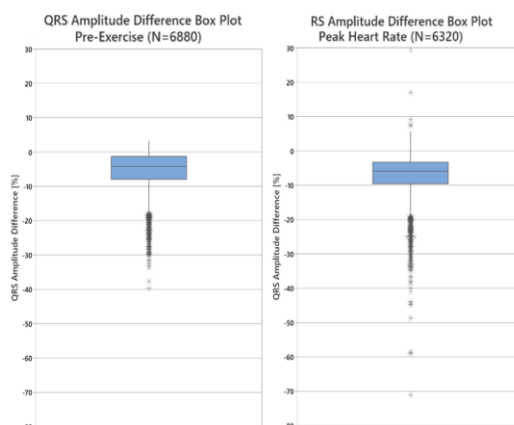


Abbildung 1: Verteilungen der QRS-Amplitudendifferenzmessungen ausgewertet in den 10-Schlag-Segmenten vor dem Training (links) und Spitzenherzfrequenz (rechts). Der Datensatz umfasst 102 echte Stressuntersuchungen vor und nach der Anwendung des Source Consistency Filter (SCF). Die Aufzeichnungen waren typischerweise 10-15 Minuten lang und bestanden aus 8 Abteilungen (I, II, V1-V6). Die Analyse konzentrierte sich auf zwei Datensegmente: Vor dem Training und Spitzenherzfrequenz. Die Vorübungsphase der Belastungsprüfung. Die maximale Herzfrequenz besteht aus 10 Schlägen, die mit der maximalen körperlichen Belastung des Patienten verbunden sind, und mit der daher als Worst-Case betrachtet.

Betroffene Produkte:

Artikelnummer	Produktbeschreibung	Chargennummer	UDI Nummer
Siehe Anlage 2	XScribe version 6 oder höher	Alle	Siehe Anlage 2
	Q-Stress version 6 oder höher		

Mögliche Gefahr:

Wenn sich die angezeigte QRS-Amplitude verändert und vom Arzt nicht erkannt wird, kann ein medizinischer Eingriff unterlassen, verzögert oder vorausgesetzt werden, der im Widerspruch zum tatsächlichen Zustand des Patienten steht.

Wenn der SCF aktiviert ist, können Diagnosen, die eine genaue Darstellung der QRS-Amplitude in der Wellenform erfordern, falsch sein, wenn sie auf der Grundlage der EKG-Wellenform erstellt werden, wie sie auf der (1) Echtzeitanzeige auf dem Bildschirm, (2) Live-EKG-Ausdrucken und (3) Abschlussberichte.

Es gibt keine klinisch signifikanten Änderungen an anderen Aspekten der Wellenform oder des angezeigten Inhalts/Abschlussberichts (z.B. dem ST-Segment, berechneten Messungen, durchschnittlichen Herzschlägen); Daher sind darauf basierende Diagnosen nicht betroffen.

Baxter identifizierte eine Kundenbeschwerde, die auf niedrige Amplitude, kleine oder falsche QRS-Messungen hinwies. Es gab keine Beschwerden im Zusammenhang mit Verletzungen oder Todesfällen.

Von Baxter zu ergreifende Massnahmen:

Aufgrund der möglichen Auswirkungen auf den QRS-Komplex entwickelt Baxter ein Software-Update für Q-Stress Version 6 oder höher und XScribe Version 6 oder höher, um das Problem zu beheben. Baxter wird Sie kontaktieren, sobald das Software-Update verfügbar ist.

Von Ihnen durchzuführende Massnahmen:

Baxter bittet Sie freundlich, die folgenden Massnahmen zu ergreifen:

- 1) Bis ein Software-Update verfügbar und installiert ist, sollten Ärzte die potentiellen Auswirkungen des SCF wie beschrieben bewerten, um zu bestimmen, ob es während des Stresstests aktiviert oder deaktiviert werden sollte.
- 2) Aktivieren oder deaktivieren Sie den SCF vor dem Stresstest durch
 - Navigieren Sie zum Menü „Modalitätseinstellungen“ (siehe Abbildung 2)

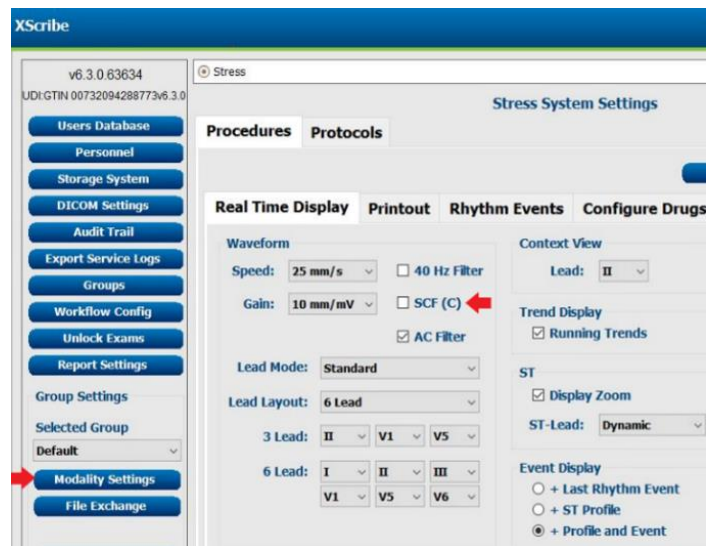


Abbildung 2 Menü „Modalitätseinstellungen“

- Es kann auch durch einen Linksklick auf eine beliebige Stelle im Echtzeit-EKG-Fenster aktiviert/deaktiviert werden, wodurch ein Fenster „Wellenformsteuerung“ geöffnet wird, in dem Benutzer die angezeigten EKG-Ableitungen, Filter (einschliesslich SCF), Anzeigeverstärkungen und Anzeigegeschwindigkeit einstellen kann (siehe Abbildung 3).

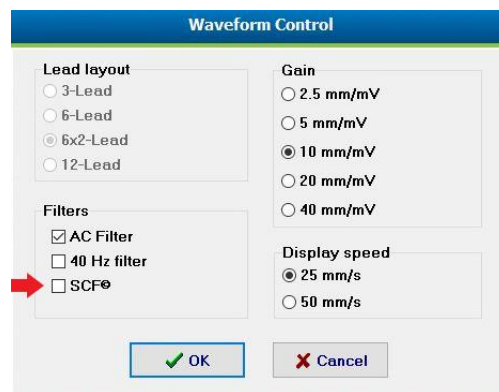


Abbildung 3 Fenster „Wellenformsteuerung“

- Wenn der Filter eingeschaltet ist, erscheint die Markierung „SCF“ am unteren rechten Rand der Echtzeit-EKG-Anzeige (siehe Abbildung 4).



Abbildung 4 Echtzeit-EKG-Anzeige

- 3) Wenn Sie dieses Produkt an andere Einrichtungen oder Abteilungen Ihrer Einrichtung verteilt haben, leiten Sie bitte eine Kopie dieser Mitteilung an diese weiter und informieren Sie sie über diese Anforderung, die potentiellen Auswirkungen des SCF während Stresstests zu bewerten
- 4) Füllen Sie das beigefügte Kundenantwortformular aus und senden Sie es an Baxter zurück, indem Sie es entweder faxen (044 908 50 40) oder einscannen und per E-Mail an (Switzerland_SHS_CQA@baxter.com) senden. Wenn Sie das Antwortformular umgehend zurücksenden, bestätigen Sie den Erhalt dieser Benachrichtigung und verhindern, dass wir Sie erneut benachrichtigen. **Dieser Schritt ist gemäss den Aufsichtsbehörden erforderlich.**
- 5) Wenn Sie dieses Produkt von einem Händler, Grosshändler, Distributor/Wiederverkäufer oder Originalgerätehersteller (OEM) sind, die betroffenen Produkte an andere Einrichtungen vertrieben hat, informieren Sie bitte ihre Kunden gemäss Ihren üblichen Verfahren über diese Korrektur.
- 6) Wenn Sie dieses Produkt von einem Händler gekauft haben, beachten Sie bitte, dass das Antwortformular von Baxter nicht anwendbar ist. Wenn Ihr Händler oder Grosshändler Ihnen ein Antwortformular zur Verfügung stellt, senden Sie es bitte gemäss den Anweisungen des Lieferanten zurück.

Die Swissmedic wurde über diese Massnahme in Kenntnis gesetzt.

Bei weiteren Fragen zu dieser Mitteilung rufen Sie uns bitte unter folgender Telefonnummer an: 044 908 50 50 oder senden Sie uns eine E-Mail an Switzerland_SHS_CQA@baxter.com.

Wir entschuldigen uns für etwaige Unannehmlichkeiten, die Ihnen und Ihren Mitarbeitern dadurch entstehen können.

Freundliche Grüsse

Baxter AG

Larissa Meyer, CQA Lead Switzerland, fvP

Anlagen:

- 1) Rückantwortformular
- 2) Liste der betroffenen Produkte

Rückantwortformular

Korrektur vom 20. Dezember 2023

Produkt: Q-stress 6 und Xscribe 6

Artikelnummer: siehe Anlage 2

Chargennummer: siehe Anlage 2

Bitte senden Sie dieses Formular ausgefüllt an die folgende Faxnummer / E-Mail Adresse:

Fax Nr. 044 908 50 40

oder

Switzerland_SHS_CQA@baxter.com

Name und Anschrift der Einrichtung:

Bitte listen Sie unten die spezifischen Produkte und Chargennummern in Ihrer Einrichtung auf*

Artikelnummer	Chargennummer

*Sie können bei Bedarf ein zusätzliches Blatt beifügen

Mit Ihrer Unterschrift bestätigen Sie, dass Sie das beigefügte Schreiben erhalten haben; führen Sie die im Schreiben beschriebenen Massnahmen nach Bedarf durch und geben Sie diese Informationen gegebenenfalls an das Personal und andere Dienste oder Einrichtungen weiter

Ausgefüllt von: _____

(bitte Name in Druckschrift)

Funktion: _____

(bitte Position/Funktion eintragen)

Telefonnummer: _____

Unterschrift: _____

Datum: ____ / ____ / ____

Liste der betroffenen Produkte

Artikelnummer	Produktbeschreibung	UDI Nummer
41000-030-50	X-SCRIBE CPI KIT WIRELESS USB 608 MHz	N/A
41000-030-51	X-SCRIBE CPI KIT WIRELESS USB 2500 MHz	N/A
41000-030-52	X-SCRIBE CPI KIT DIRECT CONNECT USB	N/A
41000-030-53	X-SCRIBE CPI KIT QUICK CONNECT USB	N/A
41000-030-60	XSCRIBE CP KIT WIRELESS	00817655020051
41000-030-62	XSCRIBE CP KIT WIRED	00817655020068
41000-030-64	XSCRIBE CP SW KIT	00817655021423
QR-6AC-AXXXX	QS6 REVIEW WKSTN KIT	00732094263961
QR-6AC-XXXXX	QS6 REVIEW WKSTN KIT	00812345027686
QS-6AA-ACAAB	QS6 AM12Q TTL SYS Z200 PCRT AHA DICOM	00732094263947
QS-6AA-CDCAA	QS6 AM12Q TTL SYS Z200 PCRT IEC XML	00732094286021
QS-6AA-DDAAB	QS6 AM12Q TTL SYS Z200 PCRT AHA DICOM	00732094263930
QS-6AA-DDCAA	QS6 AM12Q TTL SYS Z200 PCRT IEC XML	00817655021355
QS-6AA-DEAXX	QS6 AM12Q TTL SYS Z200 BCRT AHA	00732094263923
QS-6AA-XCAAA	QS6 AM12Q TTL SYS Z200 PCRT AHA XML	00732094326710
QS-6AB-BXAXA	QS6 AM12Q TTL SYS AHA XML	00732094263916
QS-6AB-CEAAA	QS6 AM12Q TTL SYS Z200 BCRT AHA XML	00732094263909
QS-6AB-CGAAB	QS6 AM12Q TTL SYS BCRT AHA DICOM	00732094263893
QS-6AC-ACAAB	QS6 AM12Q TTL KIT Z200 PCRT AHA DICOM	00732094263879
QS-6AC-AFCXA	QS6 AM12Q TTL KIT Z200 BCRT IEC XML	00732094263862
QS-6AC-AXAAA	QS6 AM12Q TTL KIT AHA XML	00732094263855
QS-6AC-CDBXA	QS6 AM12Q TTL KIT Z200 PCRT AHA XML	00732094263831
QS-6AC-CECAA	QS6 AM12Q TTL KIT Z200 BCRT IEC XML	00732094336801
QS-6AC-CFAAA	QS6 AM12Q TTL KIT Z200 BCRT AHA XML	00732094263824
QS-6AC-CFBAA	QS6 AM12Q TTL KIT Z200 BCRT AHA XML	00732094318111
QS-6AC-CFBAB	QS6 AM12Q TTL KIT Z200 BCRT AHA DICOM	00732094318029
QS-6AC-CFCAA	QS6 AM12Q TTL KIT Z200 BCRT IEC XML	00732094263817
QS-6AC-CFEAA	QS6 AM12Q TTL KIT Z200 BCRT IEC XML	00732094263800
QS-6AC-CGAAB	QS6 AM12Q TTL KIT BCRT AHA DICOM	00732094326604
QS-6AC-CXBXA	QS6 AM12Q TTL KIT AHA XML	00732094263794
QS-6AC-DDBAA	QS6 AM12Q TTL KIT Z200 PCRT AHA XML	00732094312973
QS-6AC-DDBAB	QS6 AM12Q TTL KIT Z200 PCRT AHA DICOM	00732094312966
QS-6AC-DDBAC	QS6 AM12Q TTL KIT Z200 PCRT AHA SVR	00732094312959
QS-6AC-DEAAX	QS6 AM12Q TTL KIT Z200 BCRT AHA	00732094263787
QS-6AC-DFAAA	QS6 AM12Q TTL KIT Z200 BCRT AHA XML	00732094263770
QS-6AC-DFAAX	QS6 AM12Q TTL KIT Z200 BCRT AHA	00732094263763
QS-6AC-DFBAA	QS6 AM12Q TTL KIT Z200 BCRT AHA XML	00732094302073
QS-6AC-DXAXA	QS6 AM12Q TTL KIT AHA XML	00732094263756
QS-6AC-EEBXX	QS6 AM12Q TTL KIT Z200 BCRT AHA	N/A
QS-6AC-XAAXA	QS6 AM12Q TTL KIT Z200 AHA XML	00732094263749
QS-6AC-XXAXA	QS6 AM12Q TTL KIT AHA XML	00732094355314
QS-6AC-XXAXX	QS6 AM12Q TTL KIT AHA	00732094301076
QS-6AC-XXBXA	QS6 AM12Q TTL KIT AHA XML	00732094263732
QS-6AC-XXBXX	QS6 AM12Q TTL KIT AHA	00732094346633
QS-6AD-AEAAA	QS6 AM12Q TTL SYS Z200 BCRT AHA XML	00732094288940

QS6-ATTD1	QS6 AM12Q TTL TMX SYS Z200 PCRT AHA DIC	0081234502756321
-----------	---	------------------

Artikelnummer	Produktbeschreibung	UDI Nummer
QS6-ATTDX	QS6 AM12Q TTL SYS Z200 PCRT AHA DICOM	0081234502757021
QS6-AXTDS	QS6 AM12Q TTL KIT AHA DICOM	00812345027587
QS-6BA-ACAAA	QS6 AM12Q TTL TMX SYS Z200 PCRT AHA XML	00732094263688
QS-6AD-BCAAX	QS6 AM12Q TTL SYS Z200 PCRT AHA	00732094263718
QS-6BA-ACAAB	QS6 AM12Q TTL TMX SYS Z200 PCRT AHA DIC	00732094263671
QS-6BA-ACBAA	QS6 AM12Q TTL TMX SYS Z200 PCRT AHA XML	00732094263664
QS-6BA-DDDAA	QS6 AM12Q TTL TMX SYS Z200 PCRT AHA XML	00732094263640
QS-6BB-CFAAB	QS6 AM12Q TTL TMX SYS Z200 BCRT AHA DIC	00732094263626
QS-6BC-CDCXA	QS6 AM12Q TTL TMX KIT Z200 PCRT IEC XML	00732094263596
QS-6BC-CGAAA	QS6 AM12Q TTL TMX KIT BCRT AHA XML	00732094263589
QS-6BD-ACAAB	QS6 AM12Q TTL TMX SYS Z200 PCRT AHA DIC	00812345027594
QS6-BLXC1	QS6 AM12Q TMX SYS Z200 PCRT AHA XML	0081234502760021
QS6-BLXCX	QS6 AM12Q SYS Z200 PCRT AHA XML	0081234502761721
QS6-BXXCS	QS6 AM12Q KIT AHA XML	00732094263176
QS-6CA-DFCAA	QS6 AM12QSYS Z200 BCRT IEC XML	00732094263565
QS-6CA-DGCAA	QS6 AM12QSYS BCRT IEC XML	00732094263558
QS-6CB-CXCXA	QS6 AM12QSYS IEC XML	00732094263534
QS-6CB-DFCAX	QS6 AM12QSYS Z200 BCRT IEC	00732094263527
QS-6CC-DFCAX	QS6 AM12QKIT Z200 BCRT IEC	00732094263497
QS-6CC-DFCXX	QS6 AM12Q KIT Z200 BCRT IEC	00732094308723
QS-6CC-DGCAX	QS6 AM12QKIT BCRT IEC	00732094263480
QS-6CC-DGCXX	QS6 AM12Q KIT BCRT IEC	00732094308716
QS-6CC-XAAXX	QS6 AM12QKIT Z200 AHA	00732094324075
QS-6CC-XGAAB	QS6 AM12QKIT BCRT AHA DICOM	00732094323726
QS-6CC-XXAXA	QS6 AM12QKIT AHA XML	00732094313819
QS-6CC-XXBXX	QS6 AM12QKIT AHA	00732094263466
QS-6EA-ACDAA	QS6 WAM TTL SYS Z200 PCRT AHA XML	00732094302110
QS-6EA-CDEAB	QS6 WAM TTL SYS Z200 PCRT IEC DICOM	00732094263435
QS-6EA-DFEAB	QS6 WAM TTL SYS Z200 BCRT IEC DICOM	00732094263428
QS-6EB-BXDAX	QS6 WAM TTL SYS AHA XML	00812345027693
QS-6EB-CFDAB	QS6 WAM TTL SYS Z200 BCRT AHA DICOM	00732094263411
QS-6EB-DGEAB	QS6 WAM TTL SYS BCRT IEC DICOM	00732094315660
QS-6EB-DGEAC	QS6 WAM TTL SYS BCRT IEC SVR	00732094296167
QS-6EB-DGEXX	QS6 WAM TTL SYS BCRT IEC	00732094321203
QS-6EB-DXDXX	QS6 WAM TTL SYS AHA	00732094263404
QS-6EB-DXEXB	QS6 WAM TTL SYS IEC DICOM	00732094263398
QS-6EC-ACAAB	QS6 WAM TTL KIT Z200 PCRT AHA DICOM	00732094263381
QS-6EC-ACDAB	QS6 WAM TTL KIT Z200 PCRT AHA DICOM	00732094263374
QS-6EC-ACDAX	QS6 WAM TTL KIT Z200 PCRT AHA	00732094263367
QS-6EC-AXDXA	QS6 WAM TTL KIT AHA XML	00732094263350
QS-6EC-CFEAA	QS6 WAM TTL KIT Z200 BCRT IEC XML	00732094318463
QS-6EC-DBEAB	QS6 WAM TTL KIT Z200 IEC DICOM	00732094312942
QS-6EC-DFDAA	QS6 WAM TTL KIT Z200 BCRT AHA XML	00732094263343
QS-6EC-DFEAB	QS6 WAM TTL KIT Z200 BCRT IEC DICOM	00732094312928

QS-6EC-DFEAC	QS6 WAM TTL KIT Z200 BCRT IEC SVR	00732094318074
QS-6EC-DGEAA	QS6 WAM TTL KIT BCRT IEC XML	00732094294026

Artikelnummer	Produktbeschreibung	UDI Nummer
QS-6EC-DGEAB	QS6 WAM TTL KIT BCRT IEC DICOM	00732094343786
QS-6EC-DXEXB	QS6 WAM TTL KIT IEC DICOM	00732094263336
QS-6EC-XGDAA	QS6 WAM TTL KIT BCRT AHA XML	00732094263329
QS-6EC-XXDXA	QS6 WAM TTL KIT AHA XML	00732094263312
QS-6EC-XXDXX	QS6 WAM TTL KIT AHA	00732094263305
QS-6ED-AEAAA	QS6 WAM TTL SYS Z200 BCRT AHA XML	00732094355543
QS-6ED-AEAAB	QS6 WAM TTL SYS Z200 BCRT AHA DICOM	00812345027624
QS-6ED-CBEAA	QS6 WAM TTL SYS Z200 IEC XML	00732094263299
QS-6ED-CDEAB	QS6 WAM TTL SYS Z200 PCRT IEC DICOM	00732094309645
QS-6ED-CFEAA	QS6 WAM TTL SYS Z200 BCRT IEC XML	00732094313918
QS-6ED-CFEAB	QS6 WAM TTL SYS Z200 BCRT IEC DICOM	00732094318432
QS-6ED-DFEAB	QS6 WAM TTL SYS Z200 BCRT IEC DICOM	00732094312904
QS-6ED-XEDAA	QS6 WAM TTL SYS Z200 BCRT AHA XML	00732094308495
QS-6ED-XFDAB	QS6 WAM TTL SYS Z200 BCRT AHA DICOM	00732094352672
QS-6FA-AXDAA	QS6 WAM SYS AHA XML	00732094321012
QS-6FA-CEDXA	QS6 WAM SYS Z200 BCRT AHA XML	00732094353587
QS-6FA-DEDXA	QS6 WAM SYS Z200 BCRT AHA XML	00732094301717
QS-6FB-DGEXX	QS6 WAM SYS BCRT IEC	00732094307245
QS-6FB-DXEXX	QS6 WAM SYS IEC	00732094263282
QS-6FC-AXDXA	QS6 WAM KIT AHA XML	00817655021881
QS-6FC-DDBAB	QS6 WAM KIT Z200 PCRT AHA DICOM	00732094312898
QS-6FC-DEDAX	QS6 WAM KIT Z200 BCRT AHA	00732094263275
QS-6FC-DXDAA	QS6 WAM KIT AHA XML	00732094285994
QS-6FC-DXDAX	QS6 WAM KIT AHA	00732094263251
QS-6FC-XXDAA	QS6 WAM KIT AHA XML	00732094299977
QS-6FC-XXDXA	QS6 WAM KIT AHA XML	00732094313802
QS-6GA-ACDAA	QS6 WAM TTL RS SYS Z200 PCRT AHA XML	00732094263244
QS-6GA-CDEAA	QS6 WAM TTL RS SYS Z200 PCRT IEC XML	00732094300673
QS-6GA-DDEAA	QS6 WAM TTL RS SYS Z200 PCRT IEC XML	00732094312874
QS-6GC-DXEXB	QS6 WAM TTL RS KIT IEC DICOM	00732094263237
QS-6GD-AFDAX	QS6 WAM TTL RS SYS Z200 BCRT AHA	00732094263220
QS-6GD-AGDAA	QS6 WAM TTL RS SYS BCRT AHA XML	00732094319934
QS-6GD-CAEAB	QS6 WAM TTL RS SYS Z200 IEC DICOM	00732094289169
QS-6GD-CDEAB	QS6 WAM TTL RS SYS Z200 PCRT IEC DICOM	00732094308426
QS-6GD-CFEAA	QS6 WAM TTL RS SYS Z200 BCRT IEC XML	00732094318425
QS-6GD-XGDAA	QS6 WAM TTL RS SYS BCRT AHA XML	00817655021331
QS-6IA-AEDAB	QS6 WAM RS SYS Z200 BCRT AHA DICOM	00817655020754
QS-6IA-XGDAA	QS6 WAM RS SYS BCRT AHA XML	00732094307818
QS-6IB-CGEXB	QS6 WAM RS SYS BCRT IEC DICOM	00732094263213
QS6-MLTC1	QS6 AM12Q TTL TMX SYS Z200 PCRT AHA XML	0081234502763121
QS6-MLTCX	QS6 AM12Q TTL SYS Z200 PCRT AHA XML	0081234502764821

QS6-MTTC1	QS6 AM12Q TTL TMX SYS Z200 PCRT AHA XML	0081234502765521
QS6-MTTCX	QS6 AM12Q TTL SYS Z200 PCRT AHA XML	0081234502766221
QS6-MXTCS	QS6 AM12Q TTL KIT AHA XML	0081234502767921
QS6-STTD1	QS6 PRO ENHANCED SECURITY TMILL	00732094336061
QS6-STTDX	QS6 PRO ENHANCED SECURITY	00732094336054

Artikelnummer	Produktbeschreibung	UDI Nummer
QS-6XA-DFXAA	QS6 NOACQMOD SYS Z200 BCRT NO LDSET XML	00732094263190
QS-6ID-AADAB	QS6 WAM RS SYS Z200 AHA DICOM	00732094296273
QS-6ID-AGDAA	QS6 WAM RS SYS BCRT AHA XML	00817655021935
QS-6JB-AXAXA	QS6 AM12Q RS SYS AHA XML	00732094263206
QS-6XC-BXXXX	QS6 NOACQMOD KIT NO LDSET	00732094263183
QSERV-6AA-XXXXA	QS6 SERVER KIT XML	00812345027716
QSERV-6AA-XXXXB	QS6 SERVER KIT DICOM	00812345027723
UPGXS5-B-WAM-B	UPGRADE XSCRIBE 5 WAM IEC	00732094261257
UPGXS6-B-SW-A	UPGRADE XS5 TO XS6 KIT AHA	00732094261196
XR-6AC-BXXXX	XS6/RS6 REVIEW WKSTN KIT	00732094261110
XR-6AC-DXXXX	XS6/RS6 REVIEW WKSTN KIT	00732094261103
XR-6AC-XXXXX	XS6/RS6 REVIEW WKSTN KIT	00732094261097
XSCRIBE-5AX-XXBRZ	XS6 WAM TTL KIT IEC XML BRAZIL	00817655022215
XSCRIBE-6AA-AAAAA	XS6 WAM TTL SYS Z200 AHA XML	00812345027471
QS-6XC-XXXXA	QS6 NOACQMOD KIT NO LDSET XML	00732094308457
QS-6XC-XXXXB	QS6 NOACQMOD KIT NO LDSET DICOM	00732094321029
XSCRIBE-6AA-ADBXX	XS6 WAM TTL SYS Z200 BCRT IEC	00732094261073
XSCRIBE-6AA-AFAAA	XS6 WAM TTL SYS Z200 PCRT AHA XML	00812345027488
XSCRIBE-6AA-AFAAC	XS6 WAM TTL SYS Z200 PCRT AHA SVR	00732094261066
XSCRIBE-6AA-AFAXA	XS6 WAM TTL SYS Z200 PCRT AHA XML	00732094261059
XSCRIBE-6AA-AFAXB	XS6 WAM TTL SYS Z200 PCRT AHA DICOM	00817655021997
XSCRIBE-6AA-BCAAB	XS6 WAM TTL SYS Z200 BCRT AHA DICOM	00732094299915
XSCRIBE-6AA-BXAAA	XS6 WAM TTL SYS AHA XML	00732094301489
XSCRIBE-6AA-BXAAB	XS6 WAM TTL SYS AHA DICOM	00732094261042
XSCRIBE-6AA-BXAXX	XS6 WAM TTL SYS AHA	00732094261035
XSCRIBE-6AA-CDBAA	XS6 WAM TTL SYS Z200 BCRT IEC XML	N/A
XSCRIBE-6AA-CDBAB	XS6 WAM TTL SYS Z200 BCRT IEC DICOM	00732094261028
XSCRIBE-6AA-CDBAC	XS6 WAM TTL SYS Z200 BCRT IEC SVR	00732094261011
XSCRIBE-6AA-CDBAX	XS6 WAM TTL SYS Z200 BCRT IEC	N/A
XSCRIBE-6AA-CEAAB	XS6 WAM TTL SYS BCRT AHA DICOM	00732094261004
XSCRIBE-6AA-CEBAB	XS6 WAM TTL SYS BCRT IEC DICOM	00732094260991
XSCRIBE-6AA-DEBAB	XS6 WAM TTLSYS Z200 BCRT IEC DICOM	00732094316926
XSCRIBE-6AA-DGAAA	XS6 WAM TTL SYS Z200 PCRT AHA XML	00732094260977
XSCRIBE-6AA-DGBAB	XS6 WAM TTL SYS Z200 PCRT IEC DICOM	00732094260960
XSCRIBE-6AA-DGBAX	XS6 WAM TTL SYS Z200 PCRT IEC	00732094260953
XSCRIBE-6AA-XFAAA	XS6 WAM TTL SYS Z200 PCRT AHA XML	00732094309195
XSCRIBE-6AA-XFAAB	XS6 WAM TTL SYS Z200 PCRT AHA DICOM	00732094294255
XSCRIBE-6AA-XFAXA	XS6 WAM TTL SYS Z200 PCRT AHA XML	00732094309188

XSCRIBE-6AA-XFBXA	XS6 WAM TTL SYS Z200 PCRT IEC XML	00732094309164
XSCRIBE-6AB-BEAXX	XS6 WAM TTL SYS BCRT AHA	00732094260922
XSCRIBE-6AB-BXAXA	XS6 WAM TTL SYS AHA XML	00732094260915
XSCRIBE-6AB-CDAAA	XS6 WAM TTL SYS Z200 BCRT AHA XML	00732094260908
XSCRIBE-6AB-CDAAB	XS6 WAM TTL SYS Z200 BCRT AHA DICOM	00732094260892
XSCRIBE-6AB-CDBAA	XS6 WAM TTLSYS Z200 BCRT IEC DICOM	00732094319422
XSCRIBE-6AB-CDBAB	XS6 WAM TTL SYS Z200 BCRT IEC DICOM	00732094260885
XSCRIBE-6AB-CDBAC	XS6 WAM TTL SYS Z200 BCRT IEC SVR	00732094260878

Artikelnummer	Produktbeschreibung	UDI Nummer
XSCRIBE-6AB-CEAAA	XS6 WAM TTL SYS BCRT AHA XML	00732094260861
XSCRIBE-6AB-CEBAX	XS6 WAM TTL SYS BCRT IEC	00732094260854
XSCRIBE-6AB-DDBAB	XS6 WAM TTL SYS Z200 BCRT IEC DICOM	00732094260847
XSCRIBE-6AB-DDBAX	XS6 WAM TTL SYS Z200 BCRT IEC	00732094260830
XSCRIBE-6AB-DEBAA	XS6 WAM TTL SYS BCRT IEC XML	00732094260823
XSCRIBE-6AB-DEBAB	XS6 WAM TTL SYS BCRT IEC DICOM	00732094260816
XSCRIBE-6AB-DEBAC	XS6 WAM TTL SYS BCRT IEC SVR	00732094260809
XSCRIBE-6AB-DXAAA	XS6 WAM TTL SYS AHA XML	00732094260793
XSCRIBE-6AB-DXAXA	XS6 WAM TTL SYS AHA XML	00732094260786
XSCRIBE-6AB-XCAAA	XS6 WAM TTLSYS Z200 BCRT AHA XML	00732094321081
XSCRIBE-6AC-AAAAB	XS6 WAM TTLKIT Z200 AHA DICOM	00732094339444
XSCRIBE-6AC-AAAAX	XS6 WAM TTL KIT Z200 AHA	00817655022536
XSCRIBE-6AC-AAAXX	XS6 WAM TTL KIT Z200 AHA	00817655022543
XSCRIBE-6AC-ACAAX	XS6 WAM TTL KIT Z200 BCRT AHA	00732094301120
XSCRIBE-6AC-ACAXX	XS6 WAM TTL KIT Z200 BCRT AHA	00817655021836
XSCRIBE-6AC-AEAAA	XS6 WAM TTL KIT BCRT AHA XML	00732094260779
XSCRIBE-6AC-AEAXX	XS6 WAM TTL KIT BCRT AHA	00817655022154
XSCRIBE-6AC-AXAAA	XS6 WAM TTL KIT AHA XML	00732094260762
XSCRIBE-6AC-AXAXA	XS6 WAM TTL KIT AHA XML	00812345027501
XSCRIBE-6AC-AXAXX	XS6 WAM TTL KIT AHA	00732094260755
XSCRIBE-6AC-BDAAB	XS6 WAM TTL KIT Z200 BCRT AHA DICOM	00732094260731
XSCRIBE-6AC-BEAAC	XS6 WAM TTL KIT BCRT AHA SVR	00732094260724
XSCRIBE-6AC-BEAXX	XS6 WAM TTL KIT BCRT AHA	00732094260717
XSCRIBE-6AC-BXAXA	XS6 WAM TTL KIT AHA XML	00732094260700
XSCRIBE-6AC-BXAXX	XS6 WAM TTL KIT AHA	00732094260694
XSCRIBE-6AC-CDBAA	XS6 WAM TTL KIT Z200 BCRT IEC XML	00732094260687
XSCRIBE-6AC-CDBAB	XS6 WAM TTL KIT Z200 BCRT IEC DICOM	00732094260670
XSCRIBE-6AC-CEBAX	XS6 WAM TTL KIT BCRT IEC	00732094260625
XSCRIBE-6AC-CXBAA	XS6 WAM TTL KIT IEC XML	00732094260618
XSCRIBE-6AC-CXBAB	XS6 WAM TTL KIT IEC DICOM	00732094260601
XSCRIBE-6AC-DBBXA	XS6 WAM TTL KIT Z200 IEC XML	00732094260588
XSCRIBE-6AC-DCBXA	XS6 WAM TTL KIT Z200 BCRT IEC XML	00732094260571
XSCRIBE-6AC-DDBAB	XS6 WAM TTL KIT Z200 BCRT IEC DICOM	00732094260557
XSCRIBE-6AC-DDBXB	XS6 WAM TTL KIT Z200 BCRT IEC DICOM	00732094260533
XSCRIBE-6AC-DDBAA	XS6 WAM TTL KIT Z200 BCRT IEC XML	00732094260564

XSCRIBE-6AC-DEBAA	XS6 WAM TTL KIT BCRT IEC XML	00732094260526
XSCRIBE-6AC-DEBAB	XS6 WAM TTL KIT BCRT IEC DICOM	00732094260519
XSCRIBE-6AC-DEBAC	XS6 WAM TTL KIT BCRT IEC SVR	00732094260502
XSCRIBE-6AC-DEBAX	XS6 WAM TTL KIT BCRT IEC	00732094260496
XSCRIBE-6AC-DEBXA	XS6 WAM TTL KIT BCRT IEC XML	00732094260489
XSCRIBE-6AC-DEBXX	XS6 WAM TTL KIT BCRT IEC	00732094260472
XSCRIBE-6AC-DXAAA	XS6 WAM TTL KIT AHA XML	00732094260465
XSCRIBE-6AC-DXAXX	XS6 WAM TTL KIT AHA	00732094260458
XSCRIBE-6AC-DXBAA	XS6 WAM TTL KIT IEC XML	00732094260441
XSCRIBE-6AC-DXBAB	XS6 WAM TTL KIT IEC DICOM	00732094260434
XSCRIBE-6AC-DXBAC	XS6 WAM TTL KIT IEC SVR	00732094260427

Artikelnummer	Produktbeschreibung	UDI Nummer
XSCRIBE-6AC-DXBAX	XS6 WAM TTL KIT IEC	00732094260403
XSCRIBE-6AC-DXBXA	XS6 WAM TTL KIT IEC XML	00732094260397
XSCRIBE-6AC-DXBXB	XS6 WAM TTL KIT IEC DICOM	00732094260380
XSCRIBE-6AC-DXBXD	XS6 WAM TTL KIT IEC XML CF	00732094260373
XSCRIBE-6AC-DXBXE	XS6 WAM TTL KIT IEC DICOM CF	00732094294538
XSCRIBE-6AC-DXBXX	XS6 WAM TTL KIT IEC	00732094260366
XSCRIBE-6AC-XCAXX	XS6 WAM TTLKIT Z200 BCRT AHA0	00732094308433
XSCRIBE-6AC-XDBXA	XS6 WAM TTL KIT Z200 BCRT IEC XML	00732094296587
XSCRIBE-6AC-XEBAB	XS6 WAM TTL KIT BCRT IEC DICOM	00732094260342
XSCRIBE-6AC-XXAXA	XS6 WAM TTL KIT AHA XML	00732094260335
XSCRIBE-6AC-XXAXB	XS6 WAM TTL KIT AHA DICOM	00732094301472
XSCRIBE-6AC-XXAXC	XS6 WAM TTL KIT AHA SVR	00732094301465
XSCRIBE-6AC-XXAXX	XS6 WAM TTL KIT AHA	00732094260328
XSCRIBE-6AC-XXBXX	XS6 WAM TTL KIT IEC	00732094260311
XSCRIBE-6AD-ACAAB	XS6 WAM TTLSYS Z200 BCRT AHA DICOM	00732094337518
XSCRIBE-6AD-ADAAB	XS6 WAM TTL SYS Z200 BCRT AHA DICOM	00732094260304
XSCRIBE-6AD-AEBAX	XS6 WAM TTL SYS BCRT IEC	00732094260298
XSCRIBE-6AD-AFAAA	XS6 WAM TTL SYS Z200 PCRT AHA XML	00732094260281
XSCRIBE-6AD-BCAAA	XS6 WAM TTL SYS Z200 BCRT AHA XML	00732094299892
XSCRIBE-6AD-BCAAB	XS6 WAM TTL SYS Z200 BCRT AHA DICOM	00732094299908
XSCRIBE-6AD-BEAAA	XS6 WAM TTL SYS BCRT AHA XML	00732094260274
XSCRIBE-6AD-BGAAB	XS6 WAM TTL SYS Z200 PCRT AHA DICOM	00732094302974
XSCRIBE-6AD-BXAAA	XS6 WAM TTL SYS AHA XML	00732094301458
XSCRIBE-6AD-BXAAB	XS6 WAM TTL SYS AHA DICOM	00732094296303
XSCRIBE-6AD-CDBAB	XS6 WAM TTL SYS Z200 BCRT IEC DICOM	00732094260267
XSCRIBE-6AD-CDBAC	XS6 WAM TTL SYS Z200 BCRT IEC SVR	00732094260250
XSCRIBE-6AD-DDAAA	XS6 WAM TTL SYS Z200 BCRT AHA XML	00732094260243
XSCRIBE-6AD-DBBAA	XS6 WAM TTL SYS Z200 BCRT IEC XML	00732094260236
XSCRIBE-6AD-DBBAB	XS6 WAM TTL SYS Z200 BCRT IEC DICOM	00732094260229
XSCRIBE-6AD-DEAXX	XS6 WAM TTL SYS BCRT AHA	00732094260212
XSCRIBE-6AD-DEBAA	XS6 WAM TTL SYS BCRT IEC XML	00732094260205
XSCRIBE-6AD-DEBAB	XS6 WAM TTL SYS BCRT IEC DICOM	00732094260199

XSCRIBE-6AD-XCAAA	XS6 WAM TTL SYS Z200 BCRT AHA XML	00732094309140
XSCRIBE-6AD-XCAXA	XS6 WAM TTL SYS Z200 BCRT AHA XML	00732094260182
XSCRIBE-6AD-XEAAA	XS6 WAM TTLSYS BCRT AHA XML	00732094326611
XSCRIBE-6AD-XXBAA	XS6 WAM TTL SYS IEC XML	00732094260175
XSCRIBE-6BB-DEAAB	XS6 AM12 TTL SYS BCRT AHA DICOM	00732094260151
XSCRIBE-6BB-DEBAA	XS6 AM12 TTL SYS BCRT IEC XML	00732094260144
XSCRIBE-6BB-DEBAB	XS6 AM12 TTL SYS BCRT IEC DICOM	00732094260137
XSCRIBE-6BC-AFAAB	XS6 AM12 TTL KIT Z200 PCRT AHA DICOM	00732094260120
XSCRIBE-6BC-AXAAA	XS6 AM12 TTL KIT AHA XML	00732094260113
XSCRIBE-6BC-BXAXX	XS6 AM12 TTL KIT AHA	00732094260106
XSCRIBE-6BC-CDAAA	XS6 AM12 TTL KIT Z200 BCRT AHA XML	00732094260090
XSCRIBE-6BC-DDBAA	XS6 AM12 TTL KIT Z200 BCRT IEC XML	00732094260083
XSCRIBE-6BC-DEBAA	XS6 AM12 TTL KIT BCRT IEC XML	00732094260076
XSCRIBE-6BC-DEBXX	XS6 AM12 TTL KIT BCRT IEC	00732094260069

Artikelnummer	Produktbeschreibung	UDI Nummer
XSCRIBE-6BC-DXBAA	XS6 AM12 TTL KIT IEC XML	00732094260052
XSCRIBE-6BC-DXBXX	XS6 AM12 TTL KIT IEC	00732094260045
XSCRIBE-6BC-XEBAC	XS6 AM12 TTL KIT BCRT IEC SVR	00732094260038
XSCRIBE-6BC-XXBXA	XS6 AM12 TTL KIT IEC XML	00732094260021
XSCRIBE-6BC-XXBXC	XS6 AM12 TTL KIT IEC SVR	00732094260014
XSCRIBE-6BC-XXBXD	XS6 AM12 TTL KIT IEC XML CF	00732094260007
XSCRIBE-6BC-XXBXX	XS6 AM12 TTL KIT IEC	00732094259995
XSCRIBE-6BD-CDBAB	XS6 AM12 TTL SYS Z200 BCRT IEC DICOM	00732094259971
XSCRIBE-6BD-DDAAA	XS6 AM12 TTL SYS Z200 BCRT AHA XML	00732094259957
XSCRIBE-6CA-ACAAA	XS6 WAM SYS Z200 BCRT AHA XML	00732094259940
XSCRIBE-6CA-DXAAA	XS6 WAM SYS AHA XML	00732094259926
XSCRIBE-6CB-CDBAB	XS6 WAM SYS Z200 BCRT IEC DICOM	00732094259919
XSCRIBE-6CB-DBAAA	XS6 WAM SYS Z200 AHA XML	00732094259902
XSCRIBE-6CB-DEBAA	XS6 WAM SYS BCRT IEC XML	00732094259896
XSCRIBE-6CB-DEBAB	XS6 WAM SYS BCRT IEC DICOM	00732094259889
XSCRIBE-6CB-DEBAX	XS6 WAM SYS BCRT IEC	00732094259872
XSCRIBE-6CB-DEBXX	XS6 WAM SYS BCRT IEC	00732094259865
XSCRIBE-6CB-DXBXX	XS6 WAM SYS IEC	00817655020815
XSCRIBE-6CC-AAXXX	XS6 WAM KIT Z200 NO LDSET	00732094259858
XSCRIBE-6CC-AEAAA	XS6 WAM KIT BCRT AHA XML	00732094259841
XSCRIBE-6CC-AXAAA	XS6 WAM KIT AHA XML	00732094259834
XSCRIBE-6CC-CDBAA	XS6 WAM KIT Z200 BCRT IEC XML	00732094259827
XSCRIBE-6CC-CDBAB	XS6 WAM KIT Z200 BCRT IEC DICOM	00732094259810
XSCRIBE-6CC-CDBAX	XS6 WAM KIT Z200 BCRT IEC	00732094259797
XSCRIBE-6CC-CEAAX	XS6 WAM KIT BCRT AHA	00732094259780
XSCRIBE-6CC-CEBAA	XS6 WAM KIT BCRT IEC XML	00732094259773
XSCRIBE-6CC-CEBAB	XS6 WAM KIT BCRT IEC DICOM	00732094259766
XSCRIBE-6CC-CXBAB	XS6 WAM KIT IEC DICOM	00732094259735
XSCRIBE-6CC-CXBXX	XS6 WAM KIT IEC	00732094259711

XSCRIBE-6CC-DDAAX	XS6 WAM KIT Z200 BCRT AHA	00732094259704
XSCRIBE-6CC-DDBAA	XS6 WAM KIT Z200 BCRT IEC XML	00732094259698
XSCRIBE-6CC-DDBAB	XS6 WAM KIT Z200 BCRT IEC DICOM	00732094259681
XSCRIBE-6CC-DEBAA	XS6 WAM KIT BCRT IEC XML	00732094259667
XSCRIBE-6CC-DEBAB	XS6 WAM KIT BCRT IEC DICOM	00732094259650
XSCRIBE-6CC-DEBAC	XS6 WAM KIT BCRT IEC SVR	00732094259643
XSCRIBE-6CC-DEBAD	XS6 WAM KIT BCRT IEC XML CF	00732094294552
XSCRIBE-6CC-DEBAX	XS6 WAM KIT BCRT IEC	00732094259636
XSCRIBE-6CC-DEBXA	XS6 WAM KIT BCRT IEC XML	00732094259629
XSCRIBE-6CC-DEBXX	XS6 WAM KIT BCRT IEC DICOM	00732094259612
XSCRIBE-6CC-DEBXX	XS6 WAM KIT BCRT IEC	00732094259605
XSCRIBE-6CC-DXAXX	XS6 WAM KIT AHA	00732094259599
XSCRIBE-6CC-DXBAA	XS6 WAM KIT IEC XML	00732094259582
XSCRIBE-6CC-DXBAB	XS6 WAM KIT IEC DICOM	00732094259575
XSCRIBE-6CC-DXBAX	XS6 WAM KIT IEC	00732094259551
XSCRIBE-6CC-DXBXA	XS6 WAM KIT IEC XML	00732094259544
XSCRIBE-6CC-DXBXX	XS6 WAM KIT IEC DICOM	00732094332780

Artikelnummer	Produktbeschreibung	UDI Nummer
XSCRIBE-6CC-DXBXD	XS6 WAM KIT IEC XML CF	00732094259537
XSCRIBE-6CC-DXBXE	XS6 WAM KIT IEC DICOM CF	00732094294545
XSCRIBE-6CC-DXBXX	XS6 WAM KIT IEC	00732094259520
XSCRIBE-6CC-XFAXX	XS6 WAM KIT Z200 PCRT AHA0	00732094320268
XSCRIBE-6CC-XXAXX	XS6 WAM KIT AHA	00732094259513
XSCRIBE-6CC-XXBXX	XS6 WAM KIT IEC DICOM	00732094259506
XSCRIBE-6CC-XXBXD	XS6 WAM KIT IEC XML CF	00732094259490
XSCRIBE-6CC-XXBXX	XS6 WAM KIT IEC	00732094259483
XSCRIBE-6CC-XXXXX	XS6 WAM KIT NO LDSET	00732094259476
XSCRIBE-6CD-ACAAA	XS6 WAM SYS Z200 BCRT AHA XML	00732094259469
XSCRIBE-6CD-ACAXX	XS6 WAM SYS Z200 BCRT AHA	00732094259452
XSCRIBE-6CD-CDBAB	XS6 WAM SYS Z200 BCRT IEC DICOM	00732094259445
XSCRIBE-6CD-XXBAA	XS6 WAM SYS IEC XML	00732094307825
XSCRIBE-6DB-DEBAX	XS6 AM12 SYS BCRT IEC	00732094259438
XSCRIBE-6DC-AXAAA	XS6 AM12 KIT AHA XML	00732094259421
XSCRIBE-6DC-DEBAX	XS6 AM12 KIT BCRT IEC	00732094259407
XSCRIBE-6DC-DXBAA	XS6 AM12 KIT IEC XML	00732094300383
XSCRIBE-6DC-DXBXA	XS6 AM12 KIT IEC XML	00732094259391
XSCRIBE-6DC-DXBXX	XS6 AM12 KIT IEC	00732094259384
XSCRIBE-6DC-XXBXA	XS6 AM12 KIT IEC XML	00732094259377
XSCRIBE-6DC-XXBXX	XS6 AM12 KIT IEC	00732094259360
XSCRIBE-6DD-CDBAB	XS6 AM12 SYS Z200 BCRT IEC DICOM	00732094320183
XSCRIBE-6DD-CDBAX	XS6 AM12 SYS Z200 BCRT IEC	00817655021287
XSCRIBE-6DD-DDBAB	XS6 AM12 SYS Z200 BCRT IEC DICOM	00732094320176
XSCRIBE-6EA-ACAAA	XS6 WAM TTL RS SYS Z200 BCRT AHA XML	00732094259353
XSCRIBE-6EA-DDBAB	XS6 WAM TTL RS SYS Z200 BCRT IEC DICOM	00732094259346

XSCRIBE-6EA-DEBAX	XS6 WAM TTL RS SYS BCRT IEC	00732094259339
XSCRIBE-6EA-DFBxB	XS6 WAM TTL RS SYS Z200 PCRT IEC DICOM	00732094259322
XSCRIBE-6EB-DEBAB	XS6 WAM TTL RS SYS BCRT IEC DICOM	00732094259315
XSCRIBE-6EC-DDBAB	XS6 WAM TTL RS KIT Z200 BCRT IEC DICOM	00732094259308
XSCRIBE-6EC-DEBAB	XS6 WAM TTL RS KIT BCRT IEC DICOM	00732094259292
XSCRIBE-6EC-DEBxB	XS6 WAM TTL RS KIT BCRT IEC DICOM	00732094259285
XSCRIBE-6EC-DEBXX	XS6 WAM TTL RS KIT BCRT IEC	00732094259278
XSCRIBE-6EC-DXBxB	XS6 WAM TTL RS KIT IEC DICOM	00732094259261
XSCRIBE-6EC-DXBXX	XS6 WAM TTL RS KIT IEC	00732094259254
XSCRIBE-6ED-AFAAX	XS6 WAM TTL RS SYS Z200 PCRT AHA	00732094259247
XSCRIBE-6ED-DDBAA	XS6 WAM TTL RS SYS Z200 BCRT IEC XML	00732094314168
XSCRIBE-6ED-DDBAB	XS6 WAM TTL RS SYS Z200 BCRT IEC DICOM	00732094259230
XSCRIBE-6FC-DXBXX	XS6 AM12 TTL RS KIT IEC	00732094259223
XSCRIBE-6FC-XXAXC	XS6 AM12 TTL RS KIT AHA SVR	00732094259216
XSCRIBE-6GC-DEBAA	XS6 WAM RS KIT BCRT IEC XML	00732094259186
XSCRIBE-6GC-DEBxB	XS6 WAM RS KIT BCRT IEC DICOM	00732094259179
XSCRIBE-6GC-DEBXX	XS6 WAM RS KIT BCRT IEC	00732094259162
XSCRIBE-6GC-DXBxA	XS6 WAM RS KIT IEC XML	00732094259155
XSCRIBE-6GC-DXBxB	XS6 WAM RS KIT IEC DICOM	00732094259148
XSCRIBE-6GC-DXBXX	XS6 WAM RS KIT IEC	00732094259131

Artikelnummer	Produktbeschreibung	UDI Nummer
XSCRIBE-6GC-XXBxB	XS6 WAM RS KIT IEC DICOM	00732094259100
XSCRIBE-6GD-AEAAA	XS6 WAM RS SYS BCRT AHA XML	00812345027549
XSCRIBE-6GD-CDBAB	XS6 WAM RS SYS Z200 BCRT IEC DICOM	00732094295429
XSCRIBE-6GD-DDBAX	XS6 WAM RS SYS Z200 BCRT IEC	00732094259094
XSCRIBE-6GD-XEAAA	XS6 WAM RS SYS BCRT AHA XML	00732094338485
XSCRIBE-6HB-DEBAB	XS6 AM12 RS SYS BCRT IEC DICOM	00732094259070
XSCRIBE-6XC-BXAXA	XS6 NOACQMOD KIT AHA XML	00732094259063
XSCRIBE-6XC-XXAXX	XS6 NOACQMOD KIT AHA	00732094259056
XSCRIBE-6XC-XXBXX	XS6 NOACQMOD KIT IEC0	00732094318166
XSCRIBE-6XC-XXXXA	XS6 NOACQMOD KIT NO LDSET XML	00732094259049
XSCRIBE-6XC-XXXXB	XS6 NOACQMOD KIT NO LDSET DICOM	00732094259032
XSERV-6AA-XXXXA	XS6/RS6 SERVER KIT XML	00732094257748
XSERV-6AA-XXXXB	XS6/RS6 SERVER KIT DICOM	00812345027556
XSERV-6AA-XXXXX	XS6/RS6 SERVER KIT NO CONN	00817655020730