

Opfikon, 20 Dicembre 2023

Avviso urgente di sicurezza:: FA-2023-055 Q-stress 6 e Xscribe 6


Gentile cliente,

**Descrizione
del
problema**

Baxter Healthcare Corporation sta emettendo una correzione per i sistemi per prove da sforzo cardiaco Q-Stress e XScribe (Q-Stress versione 6 o successiva e XScribe versione 6 o successiva) a causa di un potenziale cambiamento nell'ampiezza del QRS identificata nelle letture dell'elettrocardiogramma (ECG) quando il filtro (SCF: Source Consistency Filter) è abilitato. L'SCF è una funzionalità configurabile con registrazioni multicanale, che mira a ridurre il rumore e rimuovere i segnali incoerenti senza distorcere il segnale ECG associato allo stress test.

Baxter ha eseguito test interni utilizzando i dati ECG del paziente attraverso l'algoritmo SCF, che ha mostrato una variazione media dell'ampiezza QRS di una riduzione del 5,4% durante il pre-esercizio e una riduzione del 7,1% durante la frequenza cardiaca di picco (vedere Figura 1). Le variazioni massime dell'ampiezza QRS osservate sono state una riduzione del 39,7% durante il pre-esercizio e una riduzione del 71,0% durante il picco della frequenza cardiaca.

Quando l'SCF è abilitato, è possibile osservare modifiche dell'ampiezza QRS nella forma d'onda ECG visualizzata sullo schermo in tempo reale, nelle stampe ECG in tempo reale e nei rapporti finali. La visualizzazione del battito medio e tutti i calcoli (ad es. frequenza cardiaca, livello ST, pendenza ST) non si basano su dati filtrati SCF; pertanto non sono interessati. Anche gli eventi rilevati dall'algoritmo (ad esempio PVC, VRUN) non si basano su dati filtrati da SCF e pertanto non sono interessati.

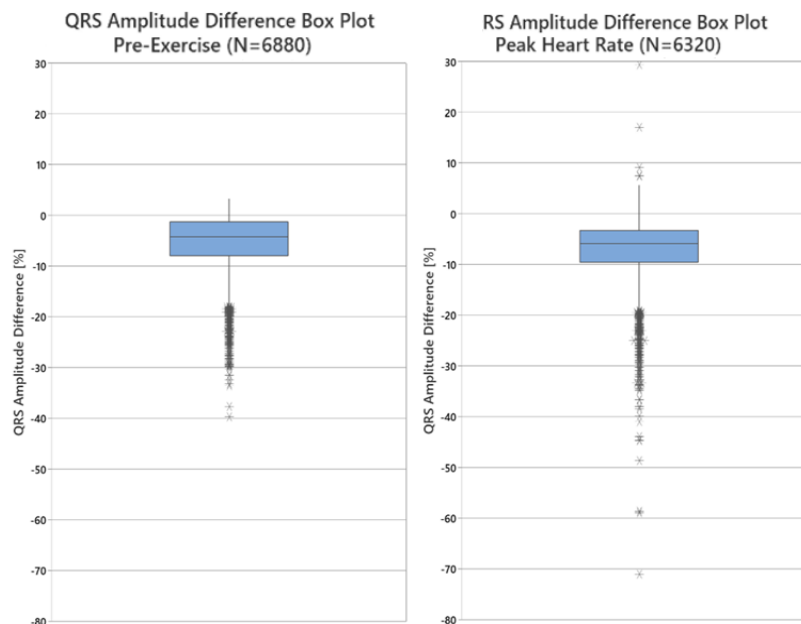**Figura 1**

Distribuzioni delle misurazioni della differenza di ampiezza QRS valutate nei segmenti di 10 battiti di Pre-esercizio (a sinistra) e Frequenza cardiaca di picco (a destra). Il set di dati include 102 esami di stress reali prima e dopo l'applicazione del Source Consistency Filter (SCF). Le registrazioni duravano tipicamente 10-15 minuti e consistevano di 8 derivazioni (I, II, V1-V6). L'analisi si è concentrata su due segmenti di dati: pre-esercizio e frequenza cardiaca di picco. Il pre-esercizio consiste nei primi 10 battiti successivi all'attivazione dell'SCF durante la fase di pre-esercizio dell'esame da sforzo. La frequenza cardiaca di picco è costituita da 10 battiti associati allo stress da sforzo massimo del paziente, quindi è considerato il caso peggiore.

**Prodotto
Impattato**

Codice prodotto	Product Description	Lot	UDI Number
Vedere allegato 2	XScribe versione 6 o successiva	tutti	Vedere allegato 2
	Q-Stress versione 6 o successiva		

**Rischi
implicati**

Se l'ampiezza QRS visualizzata viene modificata e non viene riconosciuta dal medico, l'intervento medico può essere omesso, ritardato o fornito in contraddizione con le reali condizioni del paziente.

Quando l'SCF è abilitato, le diagnosi che richiedono una rappresentazione accurata dell'ampiezza QRS nella forma d'onda potrebbero non essere corrette se effettuate in base alla forma d'onda ECG come visualizzata su (1) visualizzazione in tempo reale su schermo, (2) stampe ECG in tempo reale, e (3) relazioni finali.

Non sono presenti modifiche clinicamente significative ad altri aspetti della forma d'onda o del contenuto visualizzato/del rapporto finale (ad esempio, il segmento ST, le misurazioni calcolate, i battiti medi); pertanto, le diagnosi basate su questi non sono influenzate.

Baxter ha identificato un reclamo da parte di un cliente che indicava misurazioni QRS di bassa ampiezza, piccole o errate. Nessun reclamo è stato associato a lesioni o morte.

**Azioni da
intraprendere
da parte di
Baxter**

A causa del potenziale impatto sul complesso QRS, Baxter sta sviluppando un aggiornamento software per Q-Stress versione 6 o successiva e XScribe versione 6 o successiva per risolvere il problema. Baxter la contatterà quando l'aggiornamento del software sarà disponibile.

**Azioni da
intraprendere
da parte
dell'utente**

1. Fino a quando un aggiornamento software non sarà disponibile e non verrà installato, i medici devono valutare il potenziale impatto dell'SCF come descritto per determinare se deve essere abilitato o disabilitato durante lo stress test.
2. Abilitare o disabilitare l'SCF prima dello stress test:

- Navigazione al menu “Impostazioni Modalità ” (vedere Figura 2).

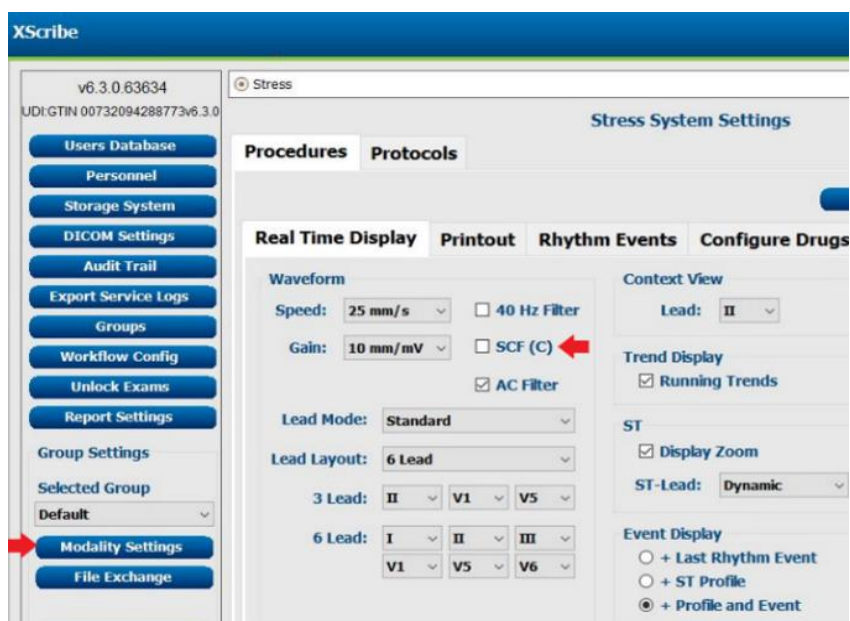


Figura 2 menu Impostazioni Modalità

- Può anche essere abilitato/disabilitato facendo clic con il pulsante sinistro del mouse in un punto qualsiasi della finestra ECG in tempo reale che apre una finestra "Controllo forma d'onda", consentendo all'utente di impostare derivazioni ECG visualizzate, filtri (incluso SCF), guadagno di visualizzazione e velocità di visualizzazione (vedere Figura 3).

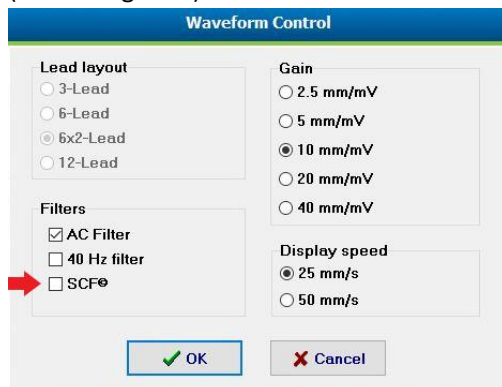


Figura 3 Finestra Controllo forma d'onda

Quando il filtro è attivo, il simbolo “SCF©” appare nel bordo inferiore destro de display ECG in tempo reale (vedere Figura 4).



Figure 4 Visualizzazione dell'ECG in tempo reale

3. Se ha distribuito questo prodotto ad altre strutture o dipartimenti della sua istituzione, inoltri loro una copia di questa comunicazione, informandoli dell'obbligo di valutare il potenziale impatto dell'SCF durante gli stress test.
4. **Se riceve questa comunicazione direttamente da Baxter, completare il modulo di risposta cliente Baxter allegato** e rispedirlo a Baxter inviandolo via fax al numero 044 908 50 40 o scansionandolo e inviandolo tramite e-mail a Switzerland_SHS_CQA@baxter.com. La restituzione tempestiva del modulo di risposta del cliente confermerà la ricezione di questa notifica e impedirà all'utente di ricevere avvisi ripetuti.
5. Se lei è un rivenditore, grossista, distributore o produttore di apparecchiature originali (OEM) che ha distribuito qualsiasi prodotto interessato ad altre strutture, la preghiamo di informare i suoi clienti di questa correzione in conformità con le procedure consuete.
6. Se il prodotto è stato acquistato da un distributore, La preghiamo di notare che il modulo di risposta cliente non è applicabile. Se il suo distributore o grossista Le invia un modulo glielo restituisca secondo le istruzioni fornite.

**Ulteriori
informazioni e
supporto**

Per eventuali necessità di chiarimenti in merito a questa comunicazione, Le chiediamo di contattare Baxter al seguente indirizzo: Switzerland_SHS_CQA@baxter.com

Il Ministero della Salute è stato informato di questa azione.

Ci scusiamo per qualsiasi inconveniente che questa comunicazione possa causare a Lei ed al Suo staff.

In fede

Baxter AG
Larissa Meyer, CQA Lead Switzerland, RP

Allegato:

- 1) modulo di risposta
- 2) Elenco dei prodotti impattati

Modulo di risposta

FA-2023-055

Comunicazione importante del Q-stress 6 e Xscribe 6

Codice: Vedere allegato 1**Lotto:** Vedere allegato 1

La preghiamo di completare il presente modulo e di trasmetterlo via fax al

numero:

044 908 50 40

oppure di inviarlo per email a:

Switzerland_SHS_CQA@baxter.com

Nome e indirizzo dell'istituzione:

Si prega di elencare di eguito i prodotti specifici e i numeri di lotto presso la vostra struttura*

numero dell'articolo	Numero di Lotto

*se necessario è possibile includere un foglio aggiuntivo

Confermiamo di aver ricevuto la lettera relativa alla suddetta comunicazione sulla sicurezza e di aver trasmesso queste informazioni alle persone e alle istituzioni interessate.

Firma: _____**Data:** ____/____/____

Annexe 2_ Elenco dei prodotti impattati

Codice prodotto	Description	UDI Number
41000-030-50	X-SCRIBE CPI KIT WIRELESS USB 608 MHz	N/A
41000-030-51	X-SCRIBE CPI KIT WIRELESS USB 2500 MHz	N/A
41000-030-52	X-SCRIBE CPI KIT DIRECT CONNECT USB	N/A
41000-030-53	X-SCRIBE CPI KIT QUICK CONNECT USB	N/A
41000-030-60	XSCRIBE CP KIT WIRELESS	00817655020051
41000-030-62	XSCRIBE CP KIT WIRED	00817655020068
41000-030-64	XSCRIBE CP SW KIT	00817655021423
QR-6AC-AXXXX	QS6 REVIEW WKSTN KIT	00732094263961
QR-6AC-XXXXX	QS6 REVIEW WKSTN KIT	00812345027686
QS-6AA-ACAAB	QS6 AM12Q TTL SYS Z200 PCRT AHA DICOM	00732094263947
QS-6AA-CDCAA	QS6 AM12Q TTL SYS Z200 PCRT IEC XML	00732094286021
QS-6AA-DDAAB	QS6 AM12Q TTL SYS Z200 PCRT AHA DICOM	00732094263930
QS-6AA-DDCAA	QS6 AM12Q TTL SYS Z200 PCRT IEC XML	00817655021355
QS-6AA-DEAXX	QS6 AM12Q TTL SYS Z200 BCRT AHA	00732094263923
QS-6AA-XCAAA	QS6 AM12Q TTL SYS Z200 PCRT AHA XML	00732094326710
QS-6AB-BXAXA	QS6 AM12Q TTL SYS AHA XML	00732094263916
QS-6AB-CEAAA	QS6 AM12Q TTL SYS Z200 BCRT AHA XML	00732094263909
QS-6AB-CGAAB	QS6 AM12Q TTL SYS BCRT AHA DICOM	00732094263893
QS-6AC-ACAAB	QS6 AM12Q TTL KIT Z200 PCRT AHA DICOM	00732094263879
QS-6AC-AFCXA	QS6 AM12Q TTL KIT Z200 BCRT IEC XML	00732094263862
QS-6AC-AXAAA	QS6 AM12Q TTL KIT AHA XML	00732094263855
QS-6AC-CDBXA	QS6 AM12Q TTL KIT Z200 PCRT AHA XML	00732094263831
QS-6AC-CECAA	QS6 AM12Q TTL KIT Z200 BCRT IEC XML	00732094336801
QS-6AC-CFAAA	QS6 AM12Q TTL KIT Z200 BCRT AHA XML	00732094263824
QS-6AC-CFBAA	QS6 AM12Q TTL KIT Z200 BCRT AHA XML	00732094318111
QS-6AC-CFBAB	QS6 AM12Q TTL KIT Z200 BCRT AHA DICOM	00732094318029
QS-6AC-CFCAA	QS6 AM12Q TTL KIT Z200 BCRT IEC XML	00732094263817
QS-6AC-CFEAA	QS6 AM12Q TTL KIT Z200 BCRT IEC XML	00732094263800
QS-6AC-CGAAB	QS6 AM12Q TTL KIT BCRT AHA DICOM	00732094326604
QS-6AC-CXBXA	QS6 AM12Q TTL KIT AHA XML	00732094263794
QS-6AC-DDBAA	QS6 AM12Q TTL KIT Z200 PCRT AHA XML	00732094312973
QS-6AC-DDBAB	QS6 AM12Q TTL KIT Z200 PCRT AHA DICOM	00732094312966
QS-6AC-DDBAC	QS6 AM12Q TTL KIT Z200 PCRT AHA SVR	00732094312959
QS-6AC-DEAAX	QS6 AM12Q TTL KIT Z200 BCRT AHA	00732094263787
QS-6AC-DFAAA	QS6 AM12Q TTL KIT Z200 BCRT AHA XML	00732094263770
QS-6AC-DFAAX	QS6 AM12Q TTL KIT Z200 BCRT AHA	00732094263763
QS-6AC-DFBAA	QS6 AM12Q TTL KIT Z200 BCRT AHA XML	00732094302073
QS-6AC-DXAXA	QS6 AM12Q TTL KIT AHA XML	00732094263756
QS-6AC-EEBXX	QS6 AM12Q TTL KIT Z200 BCRT AHA	N/A
QS-6AC-XAAXA	QS6 AM12Q TTL KIT Z200 AHA XML	00732094263749
QS-6AC-XXAXA	QS6 AM12Q TTL KIT AHA XML	00732094355314
QS-6AC-XXAXX	QS6 AM12Q TTL KIT AHA	00732094301076
QS-6AC-XXBXA	QS6 AM12Q TTL KIT AHA XML	00732094263732

QS-6AC-XXBXX	QS6 AM12Q TTL KIT AHA	00732094346633
QS-6AD-AEAAA	QS6 AM12Q TTL SYS Z200 BCRT AHA XML	00732094288940
QS6-ATTD1	QS6 AM12Q TTL TMX SYS Z200 PCRT AHA DIC	0081234502756321

Codice prodotto	Description	UDI Number
QS6-ATTDX	QS6 AM12Q TTL SYS Z200 PCRT AHA DICOM	0081234502757021
QS6-AXTDS	QS6 AM12Q TTL KIT AHA DICOM	00812345027587
QS-6BA-ACAAA	QS6 AM12Q TTL TMX SYS Z200 PCRT AHA XML	00732094263688
QS-6AD-BCAAX	QS6 AM12Q TTL SYS Z200 PCRT AHA	00732094263718
QS-6BA-ACAAB	QS6 AM12Q TTL TMX SYS Z200 PCRT AHA DIC	00732094263671
QS-6BA-ACBAA	QS6 AM12Q TTL TMX SYS Z200 PCRT AHA XML	00732094263664
QS-6BA-DDDAA	QS6 AM12Q TTL TMX SYS Z200 PCRT AHA XML	00732094263640
QS-6BB-CFAAB	QS6 AM12Q TTL TMX SYS Z200 BCRT AHA DIC	00732094263626
QS-6BC-CDCXA	QS6 AM12Q TTL TMX KIT Z200 PCRT IEC XML	00732094263596
QS-6BC-CGAAA	QS6 AM12Q TTL TMX KIT BCRT AHA XML	00732094263589
QS-6BD-ACAAB	QS6 AM12Q TTL TMX SYS Z200 PCRT AHA DIC	00812345027594
QS6-BLXC1	QS6 AM12Q TMX SYS Z200 PCRT AHA XML	0081234502760021
QS6-BLXCX	QS6 AM12Q SYS Z200 PCRT AHA XML	0081234502761721
QS6-BXXCS	QS6 AM12Q KIT AHA XML	00732094263176
QS-6CA-DFCAA	QS6 AM12QSYS Z200 BCRT IEC XML	00732094263565
QS-6CA-DGCAA	QS6 AM12QSYS BCRT IEC XML	00732094263558
QS-6CB-CXCXA	QS6 AM12QSYS IEC XML	00732094263534
QS-6CB-DFCAX	QS6 AM12QSYS Z200 BCRT IEC	00732094263527
QS-6CC-DFCAX	QS6 AM12QKIT Z200 BCRT IEC	00732094263497
QS-6CC-DFCXX	QS6 AM12Q KIT Z200 BCRT IEC	00732094308723
QS-6CC-DGCAX	QS6 AM12QKIT BCRT IEC	00732094263480
QS-6CC-DGCXX	QS6 AM12Q KIT BCRT IEC	00732094308716
QS-6CC-XAAXX	QS6 AM12QKIT Z200 AHA	00732094324075
QS-6CC-XGAAB	QS6 AM12QKIT BCRT AHA DICOM	00732094323726
QS-6CC-XXAXA	QS6 AM12QKIT AHA XML	00732094313819
QS-6CC-XXBXX	QS6 AM12QKIT AHA	00732094263466
QS-6EA-ACDAA	QS6 WAM TTL SYS Z200 PCRT AHA XML	00732094302110
QS-6EA-CDEAB	QS6 WAM TTL SYS Z200 PCRT IEC DICOM	00732094263435
QS-6EA-DFEAB	QS6 WAM TTL SYS Z200 BCRT IEC DICOM	00732094263428
QS-6EB-BXDXA	QS6 WAM TTL SYS AHA XML	00812345027693
QS-6EB-CFDAB	QS6 WAM TTL SYS Z200 BCRT AHA DICOM	00732094263411
QS-6EB-DGEAB	QS6 WAM TTL SYS BCRT IEC DICOM	00732094315660
QS-6EB-DGEAC	QS6 WAM TTL SYS BCRT IEC SVR	00732094296167
QS-6EB-DGEXX	QS6 WAM TTL SYS BCRT IEC	00732094321203
QS-6EB-DXDXX	QS6 WAM TTL SYS AHA	00732094263404
QS-6EB-DXEXB	QS6 WAM TTL SYS IEC DICOM	00732094263398
QS-6EC-ACAAB	QS6 WAM TTL KIT Z200 PCRT AHA DICOM	00732094263381
QS-6EC-ACDAB	QS6 WAM TTL KIT Z200 PCRT AHA DICOM	00732094263374
QS-6EC-ACDAX	QS6 WAM TTL KIT Z200 PCRT AHA	00732094263367
QS-6EC-AXDXA	QS6 WAM TTL KIT AHA XML	00732094263350
QS-6EC-CFEAA	QS6 WAM TTL KIT Z200 BCRT IEC XML	00732094318463
QS-6EC-DBEAB	QS6 WAM TTL KIT Z200 IEC DICOM	00732094312942
QS-6EC-DFDAA	QS6 WAM TTL KIT Z200 BCRT AHA XML	00732094263343

QS-6EC-DFEAB	QS6 WAM TTL KIT Z200 BCRT IEC DICOM	00732094312928
QS-6EC-DFEAC	QS6 WAM TTL KIT Z200 BCRT IEC SVR	00732094318074
QS-6EC-DGEAA	QS6 WAM TTL KIT BCRT IEC XML	00732094294026

Codice prodotto	Description	UDI Number
QS-6EC-DGEAB	QS6 WAM TTL KIT BCRT IEC DICOM	00732094343786
QS-6EC-DXEXB	QS6 WAM TTL KIT IEC DICOM	00732094263336
QS-6EC-XGDAA	QS6 WAM TTL KIT BCRT AHA XML	00732094263329
QS-6EC-XXDXA	QS6 WAM TTL KIT AHA XML	00732094263312
QS-6EC-XXDXX	QS6 WAM TTL KIT AHA	00732094263305
QS-6ED-AEAAA	QS6 WAM TTL SYS Z200 BCRT AHA XML	00732094355543
QS-6ED-AEAAB	QS6 WAM TTL SYS Z200 BCRT AHA DICOM	00812345027624
QS-6ED-CBEAA	QS6 WAM TTL SYS Z200 IEC XML	00732094263299
QS-6ED-CDEAB	QS6 WAM TTL SYS Z200 PCRT IEC DICOM	00732094309645
QS-6ED-CFEAA	QS6 WAM TTL SYS Z200 BCRT IEC XML	00732094313918
QS-6ED-CFEAB	QS6 WAM TTL SYS Z200 BCRT IEC DICOM	00732094318432
QS-6ED-DFEAB	QS6 WAM TTL SYS Z200 BCRT IEC DICOM	00732094312904
QS-6ED-XEDAA	QS6 WAM TTL SYS Z200 BCRT AHA XML	00732094308495
QS-6ED-XFDAB	QS6 WAM TTL SYS Z200 BCRT AHA DICOM	00732094352672
QS-6FA-AXDAA	QS6 WAM SYS AHA XML	00732094321012
QS-6FA-CEDXA	QS6 WAM SYS Z200 BCRT AHA XML	00732094353587
QS-6FA-DEDXA	QS6 WAM SYS Z200 BCRT AHA XML	00732094301717
QS-6FB-DGEXX	QS6 WAM SYS BCRT IEC	00732094307245
QS-6FB-DXEXX	QS6 WAM SYS IEC	00732094263282
QS-6FC-AXDXA	QS6 WAM KIT AHA XML	00817655021881
QS-6FC-DDBAB	QS6 WAM KIT Z200 PCRT AHA DICOM	00732094312898
QS-6FC-DEDAX	QS6 WAM KIT Z200 BCRT AHA	00732094263275
QS-6FC-DXDAA	QS6 WAM KIT AHA XML	00732094285994
QS-6FC-DXDAX	QS6 WAM KIT AHA	00732094263251
QS-6FC-XXDAA	QS6 WAM KIT AHA XML	00732094299977
QS-6FC-XXDXA	QS6 WAM KIT AHA XML	00732094313802
QS-6GA-ACDAA	QS6 WAM TTL RS SYS Z200 PCRT AHA XML	00732094263244
QS-6GA-CDEAA	QS6 WAM TTL RS SYS Z200 PCRT IEC XML	00732094300673
QS-6GA-DDEAA	QS6 WAM TTL RS SYS Z200 PCRT IEC XML	00732094312874
QS-6GC-DXEXB	QS6 WAM TTL RS KIT IEC DICOM	00732094263237
QS-6GD-AFDAX	QS6 WAM TTL RS SYS Z200 BCRT AHA	00732094263220
QS-6GD-AGDAA	QS6 WAM TTL RS SYS BCRT AHA XML	00732094319934
QS-6GD-CAEAB	QS6 WAM TTL RS SYS Z200 IEC DICOM	00732094289169
QS-6GD-CDEAB	QS6 WAM TTL RS SYS Z200 PCRT IEC DICOM	00732094308426
QS-6GD-CFEAA	QS6 WAM TTL RS SYS Z200 BCRT IEC XML	00732094318425
QS-6GD-XGDAA	QS6 WAM TTL RS SYS BCRT AHA XML	00817655021331
QS-6IA-AEDAB	QS6 WAM RS SYS Z200 BCRT AHA DICOM	00817655020754
QS-6IA-XGDAA	QS6 WAM RS SYS BCRT AHA XML	00732094307818
QS-6IB-CGEXB	QS6 WAM RS SYS BCRT IEC DICOM	00732094263213
QS6-MLTC1	QS6 AM12Q TTL TMX SYS Z200 PCRT AHA XML	0081234502763121
QS6-MLTCX	QS6 AM12Q TTL SYS Z200 PCRT AHA XML	0081234502764821
QS6-MTTC1	QS6 AM12Q TTL TMX SYS Z200 PCRT AHA XML	0081234502765521
QS6-MTTCX	QS6 AM12Q TTL SYS Z200 PCRT AHA XML	0081234502766221
QS6-MXTCS	QS6 AM12Q TTL KIT AHA XML	0081234502767921

QS6-STTD1	QS6 PRO ENHANCED SECURITY TMILL	00732094336061
QS6-STTDX	QS6 PRO ENHANCED SECURITY	00732094336054

Codice prodotto	Description	UDI Number
QS-6XA-DFXAA	QS6 NOACQMOD SYS Z200 BCRT NO LDSET XML	00732094263190
QS-6ID-AADAB	QS6 WAM RS SYS Z200 AHA DICOM	00732094296273
QS-6ID-AGDAA	QS6 WAM RS SYS BCRT AHA XML	00817655021935
QS-6JB-AXAXA	QS6 AM12Q RS SYS AHA XML	00732094263206
QS-6XC-BXXXX	QS6 NOACQMOD KIT NO LDSET	00732094263183
QSERV-6AA-XXXXA	QS6 SERVER KIT XML	00812345027716
QSERV-6AA-XXXXB	QS6 SERVER KIT DICOM	00812345027723
UPGXS5-B-WAM-B	UPGRADE XSCRIBE 5 WAM IEC	00732094261257
UPGXS6-B-SW-A	UPGRADE XS5 TO XS6 KIT AHA	00732094261196
XR-6AC-BXXXX	XS6/RS6 REVIEW WKSTN KIT	00732094261110
XR-6AC-DXXXX	XS6/RS6 REVIEW WKSTN KIT	00732094261103
XR-6AC-XXXXX	XS6/RS6 REVIEW WKSTN KIT	00732094261097
XSCRIBE-5AX-XXBRZ	XS6 WAM TTL KIT IEC XML BRAZIL	00817655022215
XSCRIBE-6AA-AAAAA	XS6 WAM TTL SYS Z200 AHA XML	00812345027471
QS-6XC-XXXXA	QS6 NOACQMOD KIT NO LDSET XML	00732094308457
QS-6XC-XXXXB	QS6 NOACQMOD KIT NO LDSET DICOM	00732094321029
XSCRIBE-6AA-ADBXX	XS6 WAM TTL SYS Z200 BCRT IEC	00732094261073
XSCRIBE-6AA-AFAAA	XS6 WAM TTL SYS Z200 PCRT AHA XML	00812345027488
XSCRIBE-6AA-AFAAC	XS6 WAM TTL SYS Z200 PCRT AHA SVR	00732094261066
XSCRIBE-6AA-AFAXA	XS6 WAM TTL SYS Z200 PCRT AHA XML	00732094261059
XSCRIBE-6AA-AFAXB	XS6 WAM TTL SYS Z200 PCRT AHA DICOM	00817655021997
XSCRIBE-6AA-BCAAB	XS6 WAM TTL SYS Z200 BCRT AHA DICOM	00732094299915
XSCRIBE-6AA-BXAAA	XS6 WAM TTL SYS AHA XML	00732094301489
XSCRIBE-6AA-BXAAB	XS6 WAM TTL SYS AHA DICOM	00732094261042
XSCRIBE-6AA-BXAXX	XS6 WAM TTL SYS AHA	00732094261035
XSCRIBE-6AA-CDBAA	XS6 WAM TTL SYS Z200 BCRT IEC XML	N/A
XSCRIBE-6AA-CDBAB	XS6 WAM TTL SYS Z200 BCRT IEC DICOM	00732094261028
XSCRIBE-6AA-CDBAC	XS6 WAM TTL SYS Z200 BCRT IEC SVR	00732094261011
XSCRIBE-6AA-CDBAX	XS6 WAM TTL SYS Z200 BCRT IEC	N/A
XSCRIBE-6AA-CEAAB	XS6 WAM TTL SYS BCRT AHA DICOM	00732094261004
XSCRIBE-6AA-CEBAB	XS6 WAM TTL SYS BCRT IEC DICOM	00732094260991
XSCRIBE-6AA-DEBAB	XS6 WAM TTLSYS Z200 BCRT IEC DICOM	00732094316926
XSCRIBE-6AA-DGAAA	XS6 WAM TTL SYS Z200 PCRT AHA XML	00732094260977
XSCRIBE-6AA-DGBAB	XS6 WAM TTL SYS Z200 PCRT IEC DICOM	00732094260960
XSCRIBE-6AA-DGBAX	XS6 WAM TTL SYS Z200 PCRT IEC	00732094260953
XSCRIBE-6AA-XFAAA	XS6 WAM TTL SYS Z200 PCRT AHA XML	00732094309195
XSCRIBE-6AA-XFAAB	XS6 WAM TTL SYS Z200 PCRT AHA DICOM	00732094294255
XSCRIBE-6AA-XFAXA	XS6 WAM TTL SYS Z200 PCRT AHA XML	00732094309188
XSCRIBE-6AA-XFBXA	XS6 WAM TTL SYS Z200 PCRT IEC XML	00732094309164
XSCRIBE-6AB-BEAXX	XS6 WAM TTL SYS BCRT AHA	00732094260922
XSCRIBE-6AB-BXAXA	XS6 WAM TTL SYS AHA XML	00732094260915
XSCRIBE-6AB-CDAAA	XS6 WAM TTL SYS Z200 BCRT AHA XML	00732094260908
XSCRIBE-6AB-CDAAB	XS6 WAM TTL SYS Z200 BCRT AHA DICOM	00732094260892
XSCRIBE-6AB-CDBAA	XS6 WAM TTLSYS Z200 BCRT IEC DICOM	00732094319422
XSCRIBE-6AB-CDBAB	XS6 WAM TTL SYS Z200 BCRT IEC DICOM	00732094260885

XSCRIBE-6AB-CDBAC	XS6 WAM TTL SYS Z200 BCRT IEC SVR	00732094260878
-------------------	-----------------------------------	----------------

Codice prodotto	Description	UDI Number
XSCRIBE-6AB-CEAAA	XS6 WAM TTL SYS BCRT AHA XML	00732094260861
XSCRIBE-6AB-CEBAX	XS6 WAM TTL SYS BCRT IEC	00732094260854
XSCRIBE-6AB-DDBAB	XS6 WAM TTL SYS Z200 BCRT IEC DICOM	00732094260847
XSCRIBE-6AB-DDBAX	XS6 WAM TTL SYS Z200 BCRT IEC	00732094260830
XSCRIBE-6AB-DEBAA	XS6 WAM TTL SYS BCRT IEC XML	00732094260823
XSCRIBE-6AB-DEBAB	XS6 WAM TTL SYS BCRT IEC DICOM	00732094260816
XSCRIBE-6AB-DEBAC	XS6 WAM TTL SYS BCRT IEC SVR	00732094260809
XSCRIBE-6AB-DXAAA	XS6 WAM TTL SYS AHA XML	00732094260793
XSCRIBE-6AB-DXAXA	XS6 WAM TTL SYS AHA XML	00732094260786
XSCRIBE-6AB-XCAAA	XS6 WAM TTLSYS Z200 BCRT AHA XML	00732094321081
XSCRIBE-6AC-AAAAB	XS6 WAM TTLKIT Z200 AHA DICOM	00732094339444
XSCRIBE-6AC-AAAAX	XS6 WAM TTL KIT Z200 AHA	00817655022536
XSCRIBE-6AC-AAAXX	XS6 WAM TTL KIT Z200 AHA	00817655022543
XSCRIBE-6AC-ACAAX	XS6 WAM TTL KIT Z200 BCRT AHA	00732094301120
XSCRIBE-6AC-ACAXX	XS6 WAM TTL KIT Z200 BCRT AHA	00817655021836
XSCRIBE-6AC-AEAAA	XS6 WAM TTL KIT BCRT AHA XML	00732094260779
XSCRIBE-6AC-AEAXX	XS6 WAM TTL KIT BCRT AHA	00817655022154
XSCRIBE-6AC-AXAAA	XS6 WAM TTL KIT AHA XML	00732094260762
XSCRIBE-6AC-AXAXA	XS6 WAM TTL KIT AHA XML	00812345027501
XSCRIBE-6AC-AXAXX	XS6 WAM TTL KIT AHA	00732094260755
XSCRIBE-6AC-BDAAB	XS6 WAM TTL KIT Z200 BCRT AHA DICOM	00732094260731
XSCRIBE-6AC-BEAAC	XS6 WAM TTL KIT BCRT AHA SVR	00732094260724
XSCRIBE-6AC-BEAXX	XS6 WAM TTL KIT BCRT AHA	00732094260717
XSCRIBE-6AC-BXAXA	XS6 WAM TTL KIT AHA XML	00732094260700
XSCRIBE-6AC-BXAXX	XS6 WAM TTL KIT AHA	00732094260694
XSCRIBE-6AC-CDBAA	XS6 WAM TTL KIT Z200 BCRT IEC XML	00732094260687
XSCRIBE-6AC-CDBAB	XS6 WAM TTL KIT Z200 BCRT IEC DICOM	00732094260670
XSCRIBE-6AC-CEBAX	XS6 WAM TTL KIT BCRT IEC	00732094260625
XSCRIBE-6AC-CXBAA	XS6 WAM TTL KIT IEC XML	00732094260618
XSCRIBE-6AC-CXBAB	XS6 WAM TTL KIT IEC DICOM	00732094260601
XSCRIBE-6AC-DBBXA	XS6 WAM TTL KIT Z200 IEC XML	00732094260588
XSCRIBE-6AC-DCBXA	XS6 WAM TTL KIT Z200 BCRT IEC XML	00732094260571
XSCRIBE-6AC-DDBAB	XS6 WAM TTL KIT Z200 BCRT IEC DICOM	00732094260557
XSCRIBE-6AC-DDBxB	XS6 WAM TTL KIT Z200 BCRT IEC DICOM	00732094260533
XSCRIBE-6AC-DDBAA	XS6 WAM TTL KIT Z200 BCRT IEC XML	00732094260564
XSCRIBE-6AC-DEBAA	XS6 WAM TTL KIT BCRT IEC XML	00732094260526
XSCRIBE-6AC-DEBAB	XS6 WAM TTL KIT BCRT IEC DICOM	00732094260519
XSCRIBE-6AC-DEBAC	XS6 WAM TTL KIT BCRT IEC SVR	00732094260502
XSCRIBE-6AC-DEBAX	XS6 WAM TTL KIT BCRT IEC	00732094260496
XSCRIBE-6AC-DEBXA	XS6 WAM TTL KIT BCRT IEC XML	00732094260489
XSCRIBE-6AC-DEBXX	XS6 WAM TTL KIT BCRT IEC	00732094260472
XSCRIBE-6AC-DXAAA	XS6 WAM TTL KIT AHA XML	00732094260465
XSCRIBE-6AC-DXAXX	XS6 WAM TTL KIT AHA	00732094260458
XSCRIBE-6AC-DXBAA	XS6 WAM TTL KIT IEC XML	00732094260441
XSCRIBE-6AC-DXBAB	XS6 WAM TTL KIT IEC DICOM	00732094260434
XSCRIBE-6AC-DXBAC	XS6 WAM TTL KIT IEC SVR	00732094260427

Codice prodotto	Description	UDI Number
XSCRIBE-6AC-DXBAX	XS6 WAM TTL KIT IEC	00732094260403
XSCRIBE-6AC-DXBXA	XS6 WAM TTL KIT IEC XML	00732094260397
XSCRIBE-6AC-DXBXB	XS6 WAM TTL KIT IEC DICOM	00732094260380
XSCRIBE-6AC-DXBXD	XS6 WAM TTL KIT IEC XML CF	00732094260373
XSCRIBE-6AC-DXBXE	XS6 WAM TTL KIT IEC DICOM CF	00732094294538
XSCRIBE-6AC-DXBXX	XS6 WAM TTL KIT IEC	00732094260366
XSCRIBE-6AC-XCAXX	XS6 WAM TTLKIT Z200 BCRT AHA0	00732094308433
XSCRIBE-6AC-XDBXA	XS6 WAM TTL KIT Z200 BCRT IEC XML	00732094296587
XSCRIBE-6AC-XEBAB	XS6 WAM TTL KIT BCRT IEC DICOM	00732094260342
XSCRIBE-6AC-XXAXA	XS6 WAM TTL KIT AHA XML	00732094260335
XSCRIBE-6AC-XXAXB	XS6 WAM TTL KIT AHA DICOM	00732094301472
XSCRIBE-6AC-XXAXC	XS6 WAM TTL KIT AHA SVR	00732094301465
XSCRIBE-6AC-XXAXX	XS6 WAM TTL KIT AHA	00732094260328
XSCRIBE-6AC-XXBXX	XS6 WAM TTL KIT IEC	00732094260311
XSCRIBE-6AD-ACAAB	XS6 WAM TTLSYS Z200 BCRT AHA DICOM	00732094337518
XSCRIBE-6AD-ADAAB	XS6 WAM TTL SYS Z200 BCRT AHA DICOM	00732094260304
XSCRIBE-6AD-AEBAX	XS6 WAM TTL SYS BCRT IEC	00732094260298
XSCRIBE-6AD-AFAAA	XS6 WAM TTL SYS Z200 PCRT AHA XML	00732094260281
XSCRIBE-6AD-BCAAA	XS6 WAM TTL SYS Z200 BCRT AHA XML	00732094299892
XSCRIBE-6AD-BCAAB	XS6 WAM TTL SYS Z200 BCRT AHA DICOM	00732094299908
XSCRIBE-6AD-BEAAA	XS6 WAM TTL SYS BCRT AHA XML	00732094260274
XSCRIBE-6AD-BGAAB	XS6 WAM TTL SYS Z200 PCRT AHA DICOM	00732094302974
XSCRIBE-6AD-BXAAA	XS6 WAM TTL SYS AHA XML	00732094301458
XSCRIBE-6AD-BXAAB	XS6 WAM TTL SYS AHA DICOM	00732094296303
XSCRIBE-6AD-CDBAB	XS6 WAM TTL SYS Z200 BCRT IEC DICOM	00732094260267
XSCRIBE-6AD-CDBAC	XS6 WAM TTL SYS Z200 BCRT IEC SVR	00732094260250
XSCRIBE-6AD-DDAAA	XS6 WAM TTL SYS Z200 BCRT AHA XML	00732094260243
XSCRIBE-6AD-DDBAA	XS6 WAM TTL SYS Z200 BCRT IEC XML	00732094260236
XSCRIBE-6AD-DDBAB	XS6 WAM TTL SYS Z200 BCRT IEC DICOM	00732094260229
XSCRIBE-6AD-DEAXX	XS6 WAM TTL SYS BCRT AHA	00732094260212
XSCRIBE-6AD-DEBAA	XS6 WAM TTL SYS BCRT IEC XML	00732094260205
XSCRIBE-6AD-DEBAB	XS6 WAM TTL SYS BCRT IEC DICOM	00732094260199
XSCRIBE-6AD-XCAAA	XS6 WAM TTL SYS Z200 BCRT AHA XML	00732094309140
XSCRIBE-6AD-XCAXA	XS6 WAM TTL SYS Z200 BCRT AHA XML	00732094260182
XSCRIBE-6AD-XEAAA	XS6 WAM TTLSYS BCRT AHA XML	00732094326611
XSCRIBE-6AD-XXBAA	XS6 WAM TTL SYS IEC XML	00732094260175
XSCRIBE-6BB-DEAAB	XS6 AM12 TTL SYS BCRT AHA DICOM	00732094260151
XSCRIBE-6BB-DEBAA	XS6 AM12 TTL SYS BCRT IEC XML	00732094260144
XSCRIBE-6BB-DEBAB	XS6 AM12 TTL SYS BCRT IEC DICOM	00732094260137
XSCRIBE-6BC-AFAAB	XS6 AM12 TTL KIT Z200 PCRT AHA DICOM	00732094260120
XSCRIBE-6BC-AXAAA	XS6 AM12 TTL KIT AHA XML	00732094260113
XSCRIBE-6BC-BXAXX	XS6 AM12 TTL KIT AHA	00732094260106
XSCRIBE-6BC-CDAAA	XS6 AM12 TTL KIT Z200 BCRT AHA XML	00732094260090
XSCRIBE-6BC-DEBAA	XS6 AM12 TTL KIT Z200 BCRT IEC XML	00732094260083
XSCRIBE-6BC-DEBAA	XS6 AM12 TTL KIT BCRT IEC XML	00732094260076
XSCRIBE-6BC-DEBXX	XS6 AM12 TTL KIT BCRT IEC	00732094260069

Codice prodotto	Description	UDI Number
XSCRIBE-6BC-DXBAA	XS6 AM12 TTL KIT IEC XML	00732094260052
XSCRIBE-6BC-DXBXX	XS6 AM12 TTL KIT IEC	00732094260045
XSCRIBE-6BC-XEBAC	XS6 AM12 TTL KIT BCRT IEC SVR	00732094260038
XSCRIBE-6BC-XXBXA	XS6 AM12 TTL KIT IEC XML	00732094260021
XSCRIBE-6BC-XXBXC	XS6 AM12 TTL KIT IEC SVR	00732094260014
XSCRIBE-6BC-XXBXD	XS6 AM12 TTL KIT IEC XML CF	00732094260007
XSCRIBE-6BC-XXBXX	XS6 AM12 TTL KIT IEC	00732094259995
XSCRIBE-6BD-CDBAB	XS6 AM12 TTL SYS Z200 BCRT IEC DICOM	00732094259971
XSCRIBE-6BD-DDAAA	XS6 AM12 TTL SYS Z200 BCRT AHA XML	00732094259957
XSCRIBE-6CA-ACAAA	XS6 WAM SYS Z200 BCRT AHA XML	00732094259940
XSCRIBE-6CA-DXAAA	XS6 WAM SYS AHA XML	00732094259926
XSCRIBE-6CB-CDBAB	XS6 WAM SYS Z200 BCRT IEC DICOM	00732094259919
XSCRIBE-6CB-DBAAA	XS6 WAM SYS Z200 AHA XML	00732094259902
XSCRIBE-6CB-DEBAA	XS6 WAM SYS BCRT IEC XML	00732094259896
XSCRIBE-6CB-DEBAB	XS6 WAM SYS BCRT IEC DICOM	00732094259889
XSCRIBE-6CB-DEBAX	XS6 WAM SYS BCRT IEC	00732094259872
XSCRIBE-6CB-DEBXX	XS6 WAM SYS BCRT IEC	00732094259865
XSCRIBE-6CB-DXBXX	XS6 WAM SYS IEC	00817655020815
XSCRIBE-6CC-AAXXX	XS6 WAM KIT Z200 NO LDSET	00732094259858
XSCRIBE-6CC-AEAAA	XS6 WAM KIT BCRT AHA XML	00732094259841
XSCRIBE-6CC-AXAAA	XS6 WAM KIT AHA XML	00732094259834
XSCRIBE-6CC-CDBAA	XS6 WAM KIT Z200 BCRT IEC XML	00732094259827
XSCRIBE-6CC-CDBAB	XS6 WAM KIT Z200 BCRT IEC DICOM	00732094259810
XSCRIBE-6CC-CDBAX	XS6 WAM KIT Z200 BCRT IEC	00732094259797
XSCRIBE-6CC-CEAAX	XS6 WAM KIT BCRT AHA	00732094259780
XSCRIBE-6CC-CEBAA	XS6 WAM KIT BCRT IEC XML	00732094259773
XSCRIBE-6CC-CEBAB	XS6 WAM KIT BCRT IEC DICOM	00732094259766
XSCRIBE-6CC-CXBAB	XS6 WAM KIT IEC DICOM	00732094259735
XSCRIBE-6CC-CXBXX	XS6 WAM KIT IEC	00732094259711
XSCRIBE-6CC-DDAAX	XS6 WAM KIT Z200 BCRT AHA	00732094259704
XSCRIBE-6CC-DDBAA	XS6 WAM KIT Z200 BCRT IEC XML	00732094259698
XSCRIBE-6CC-DDBAB	XS6 WAM KIT Z200 BCRT IEC DICOM	00732094259681
XSCRIBE-6CC-DEBAA	XS6 WAM KIT BCRT IEC XML	00732094259667
XSCRIBE-6CC-DEBAB	XS6 WAM KIT BCRT IEC DICOM	00732094259650
XSCRIBE-6CC-DEBAC	XS6 WAM KIT BCRT IEC SVR	00732094259643
XSCRIBE-6CC-DEBAD	XS6 WAM KIT BCRT IEC XML CF	00732094294552
XSCRIBE-6CC-DEBAX	XS6 WAM KIT BCRT IEC	00732094259636
XSCRIBE-6CC-DEBXA	XS6 WAM KIT BCRT IEC XML	00732094259629
XSCRIBE-6CC-DEBXB	XS6 WAM KIT BCRT IEC DICOM	00732094259612
XSCRIBE-6CC-DEBXX	XS6 WAM KIT BCRT IEC	00732094259605
XSCRIBE-6CC-DXAXX	XS6 WAM KIT AHA	00732094259599
XSCRIBE-6CC-DXBAA	XS6 WAM KIT IEC XML	00732094259582
XSCRIBE-6CC-DXBAB	XS6 WAM KIT IEC DICOM	00732094259575
XSCRIBE-6CC-DXBAX	XS6 WAM KIT IEC	00732094259551
XSCRIBE-6CC-DXBXA	XS6 WAM KIT IEC XML	00732094259544
XSCRIBE-6CC-DXBXB	XS6 WAM KIT IEC DICOM	00732094332780

Code produit	Description	UDI Number
XSCRIBE-6CC-DXBXD	XS6 WAM KIT IEC XML CF	00732094259537
XSCRIBE-6CC-DXBXE	XS6 WAM KIT IEC DICOM CF	00732094294545
XSCRIBE-6CC-DXBXX	XS6 WAM KIT IEC	00732094259520
XSCRIBE-6CC-XFAXX	XS6 WAM KIT Z200 PCRT AHA0	00732094320268
XSCRIBE-6CC-XXAXX	XS6 WAM KIT AHA	00732094259513
XSCRIBE-6CC-XXBXB	XS6 WAM KIT IEC DICOM	00732094259506
XSCRIBE-6CC-XXBXD	XS6 WAM KIT IEC XML CF	00732094259490
XSCRIBE-6CC-XXBXX	XS6 WAM KIT IEC	00732094259483
XSCRIBE-6CC-XXXXX	XS6 WAM KIT NO LDSET	00732094259476
XSCRIBE-6CD-ACAAA	XS6 WAM SYS Z200 BCRT AHA XML	00732094259469
XSCRIBE-6CD-ACAXX	XS6 WAM SYS Z200 BCRT AHA	00732094259452
XSCRIBE-6CD-CDBAB	XS6 WAM SYS Z200 BCRT IEC DICOM	00732094259445
XSCRIBE-6CD-XXBAA	XS6 WAM SYS IEC XML	00732094307825
XSCRIBE-6DB-DEBAX	XS6 AM12 SYS BCRT IEC	00732094259438
XSCRIBE-6DC-AXAAA	XS6 AM12 KIT AHA XML	00732094259421
XSCRIBE-6DC-DEBAX	XS6 AM12 KIT BCRT IEC	00732094259407
XSCRIBE-6DC-DXBAA	XS6 AM12 KIT IEC XML	00732094300383
XSCRIBE-6DC-DXBXA	XS6 AM12 KIT IEC XML	00732094259391
XSCRIBE-6DC-DXBXX	XS6 AM12 KIT IEC	00732094259384
XSCRIBE-6DC-XXBXA	XS6 AM12 KIT IEC XML	00732094259377
XSCRIBE-6DC-XXBXX	XS6 AM12 KIT IEC	00732094259360
XSCRIBE-6DD-CDBAB	XS6 AM12 SYS Z200 BCRT IEC DICOM	00732094320183
XSCRIBE-6DD-CDBAX	XS6 AM12 SYS Z200 BCRT IEC	00817655021287
XSCRIBE-6DD-DEBAB	XS6 AM12 SYS Z200 BCRT IEC DICOM	00732094320176
XSCRIBE-6EA-ACAAA	XS6 WAM TTL RS SYS Z200 BCRT AHA XML	00732094259353
XSCRIBE-6EA-DEBAB	XS6 WAM TTL RS SYS Z200 BCRT IEC DICOM	00732094259346
XSCRIBE-6EA-DEBAX	XS6 WAM TTL RS SYS BCRT IEC	00732094259339
XSCRIBE-6EA-DFBXX	XS6 WAM TTL RS SYS Z200 PCRT IEC DICOM	00732094259322
XSCRIBE-6EB-DEBAB	XS6 WAM TTL RS SYS BCRT IEC DICOM	00732094259315
XSCRIBE-6EC-DEBAB	XS6 WAM TTL RS KIT Z200 BCRT IEC DICOM	00732094259308
XSCRIBE-6EC-DEBAX	XS6 WAM TTL RS KIT BCRT IEC DICOM	00732094259292
XSCRIBE-6EC-DEBXX	XS6 WAM TTL RS KIT BCRT IEC	00732094259278
XSCRIBE-6EC-DXBXX	XS6 WAM TTL RS KIT IEC	00732094259261
XSCRIBE-6ED-AFAAX	XS6 WAM TTL RS SYS Z200 PCRT AHA	00732094259247
XSCRIBE-6ED-DEBAA	XS6 WAM TTL RS SYS Z200 BCRT IEC XML	00732094314168
XSCRIBE-6ED-DEBAB	XS6 WAM TTL RS SYS Z200 BCRT IEC DICOM	00732094259230
XSCRIBE-6FC-DXBXX	XS6 AM12 TTL RS KIT IEC	00732094259223
XSCRIBE-6FC-XXAXC	XS6 AM12 TTL RS KIT AHA SVR	00732094259216
XSCRIBE-6GC-DEBAA	XS6 WAM RS KIT BCRT IEC XML	00732094259186
XSCRIBE-6GC-DEBXX	XS6 WAM RS KIT BCRT IEC	00732094259179
XSCRIBE-6GC-DXBXA	XS6 WAM RS KIT IEC XML	00732094259162
XSCRIBE-6GC-DXBXX	XS6 WAM RS KIT IEC	00732094259155
XSCRIBE-6GC-DXBXB	XS6 WAM RS KIT IEC DICOM	00732094259148
XSCRIBE-6GC-DXBXX	XS6 WAM RS KIT IEC	00732094259131

Codice prodotto	Description	Lot UDI Number
XSCRIBE-6GC-XXBXB	XS6 WAM RS KIT IEC DICOM	00732094259100
XSCRIBE-6GD-AEAAA	XS6 WAM RS SYS BCRT AHA XML	00812345027549
XSCRIBE-6GD-CDBAB	XS6 WAM RS SYS Z200 BCRT IEC DICOM	00732094295429
XSCRIBE-6GD-DDBAX	XS6 WAM RS SYS Z200 BCRT IEC	00732094259094
XSCRIBE-6GD-XEAAA	XS6 WAM RS SYS BCRT AHA XML	00732094338485
XSCRIBE-6HB-DEBAB	XS6 AM12 RS SYS BCRT IEC DICOM	00732094259070
XSCRIBE-6XC-BXAXA	XS6 NOACQMOD KIT AHA XML	00732094259063
XSCRIBE-6XC-XXAXX	XS6 NOACQMOD KIT AHA	00732094259056
XSCRIBE-6XC-XXBXX	XS6 NOACQMOD KIT IEC0	00732094318166
XSCRIBE-6XC-XXXXA	XS6 NOACQMOD KIT NO LDSET XML	00732094259049
XSCRIBE-6XC-XXXXB	XS6 NOACQMOD KIT NO LDSET DICOM	00732094259032
XSERV-6AA-XXXXA	XS6/RS6 SERVER KIT XML	00732094257748
XSERV-6AA-XXXXB	XS6/RS6 SERVER KIT DICOM	00812345027556
XSERV-6AA-XXXXX	XS6/RS6 SERVER KIT NO CONN	00817655020730